

Canadian Bioinformatics Workshops

www.bioinformatics.ca

bioinformaticsdotca.github.io

Supported by



This page is available in the following languages:

Afrikaans Български Català Dansk Deutsch Ελληνικά English English (CA) English (GB) English (US) Esperanto
Castellano Castellano (AR) Español (CL) Castellano (CO) Español (Ecuador) Castellano (MX) Castellano (PE)
Euskara Suomeksi français français (CA) Galego עברית hrvatski Magyar Italiano 日本語 한국어 Macedonian Melayu
Nederlands Norsk Sesotho sa Leboa póski Português română slovenski jezik čeština srpski (latinica) Sotho svenska
中文 漢語 (台灣) isiZulu

 creative commons

Attribution-Share Alike 2.5 Canada

You are free:

 to Share — to copy, distribute and transmit the work

 to Remix — to adapt the work





Under the following conditions:

 **Attribution.** You must attribute the work in the manner specified by the author or licensor (but not in any way that suggests that they endorse you or your use of the work).

 **Share Alike.** If you alter, transform, or build upon this work, you may distribute the resulting work only under the same or similar licence to this one.

- For any reuse or distribution, you must make clear to others the licence terms of this work.
- Any of the above conditions can be waived if you get permission from the copyright holder.
- The author's moral rights are retained in this licence.

[Disclaimer](#)

Your fair dealing and other rights are in no way affected by the above.
This is a human-readable summary of the Legal Code (the full licence) available in the following languages:
[English](#) [French](#)

Learn how to distribute your work using this licence



Introduction to IGV: The Integrative Genomics Viewer

Kelsy Cotto, Obi Griffith, Malachi Griffith, Saad Khan, Allegra Petti, Huiming Xia

Informatics for RNA-seq Analysis

June 17-19, 2020



Visualization Tools in Genomics

- there are **over 40 different genome browsers**, which to use?
- depends on
 - task at hand
 - kind and size of data
 - data privacy

HT-seq Genome Browsers



Integrative
Genome
Viewer



UCSC
Genome Browser
Cancer Genome Browser



Trackster
(part of Galaxy)




Savant
Genome
Browser

- task at hand : visualizing HT-seq reads, especially good for inspecting variants
- kind and size of data : large BAM files, stored locally or remotely
- data privacy : run on the desktop, can keep all data private
- UCSC Genome Browser has been retro-fitted to display BAM files
- Trackster is a genome browser that can perform visual analytics on small windows of the genome, deploy full analysis with Galaxy

Integrative Genomics Viewer (IGV)

Desktop application for the interactive visual exploration of integrated genomic datasets

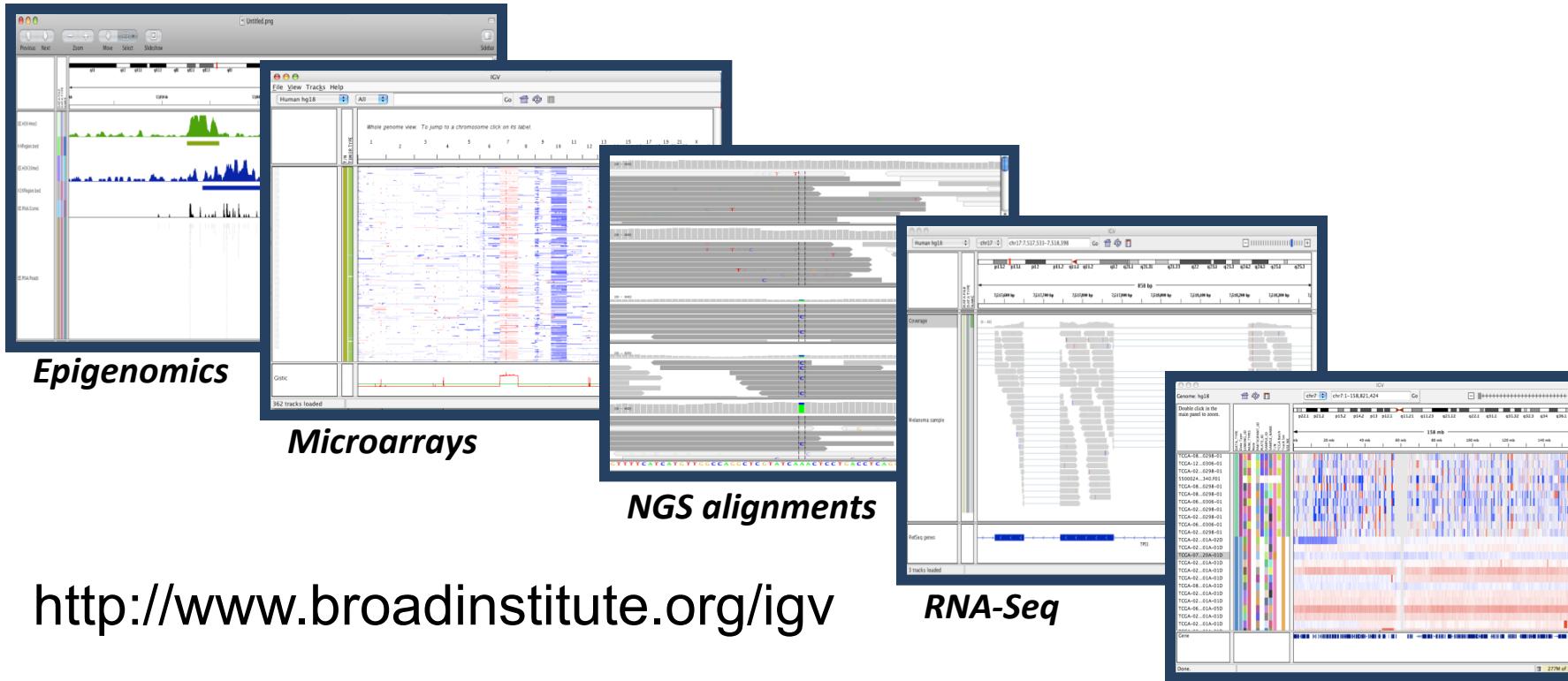

<http://www.broadinstitute.org/igv>


>85,000 registrations (2014)

Features

With IGV you can...

- Explore large genomic datasets with an intuitive, easy-to-use interface.
- Integrate multiple data types with clinical and other sample information.
- View data from multiple sources:
 - local, remote, and “cloud-based”.
- Automation of specific tasks using command-line interface

IGV data sources

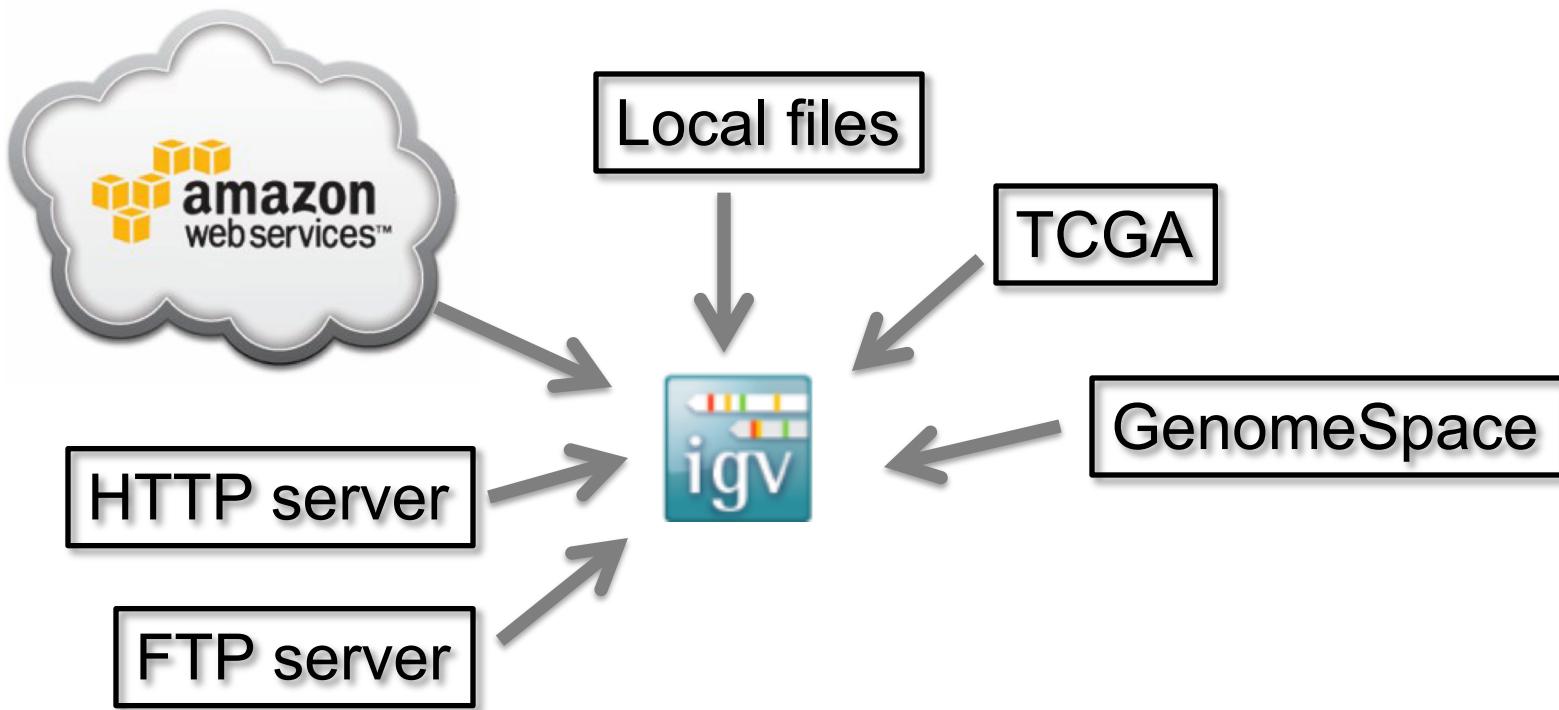


- View **local** files without uploading.
- View **remote** files without downloading the whole dataset.


Using IGV: the basics

- Launch IGV
- Select a reference genome
- Load data
- Navigate through the data
 - WGS data
 - SNVs
 - structural variations


Launch IGV

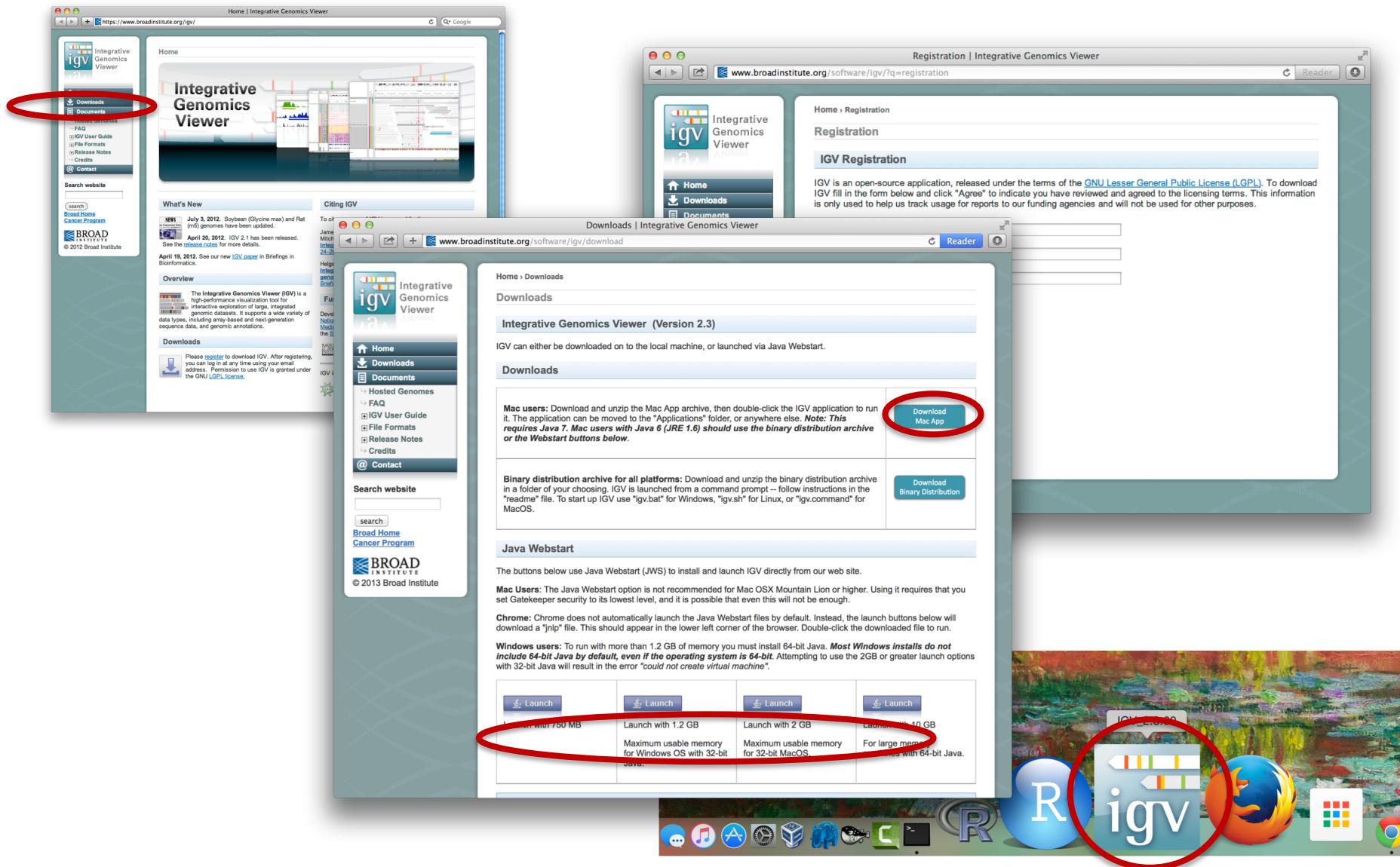


Launch IGV

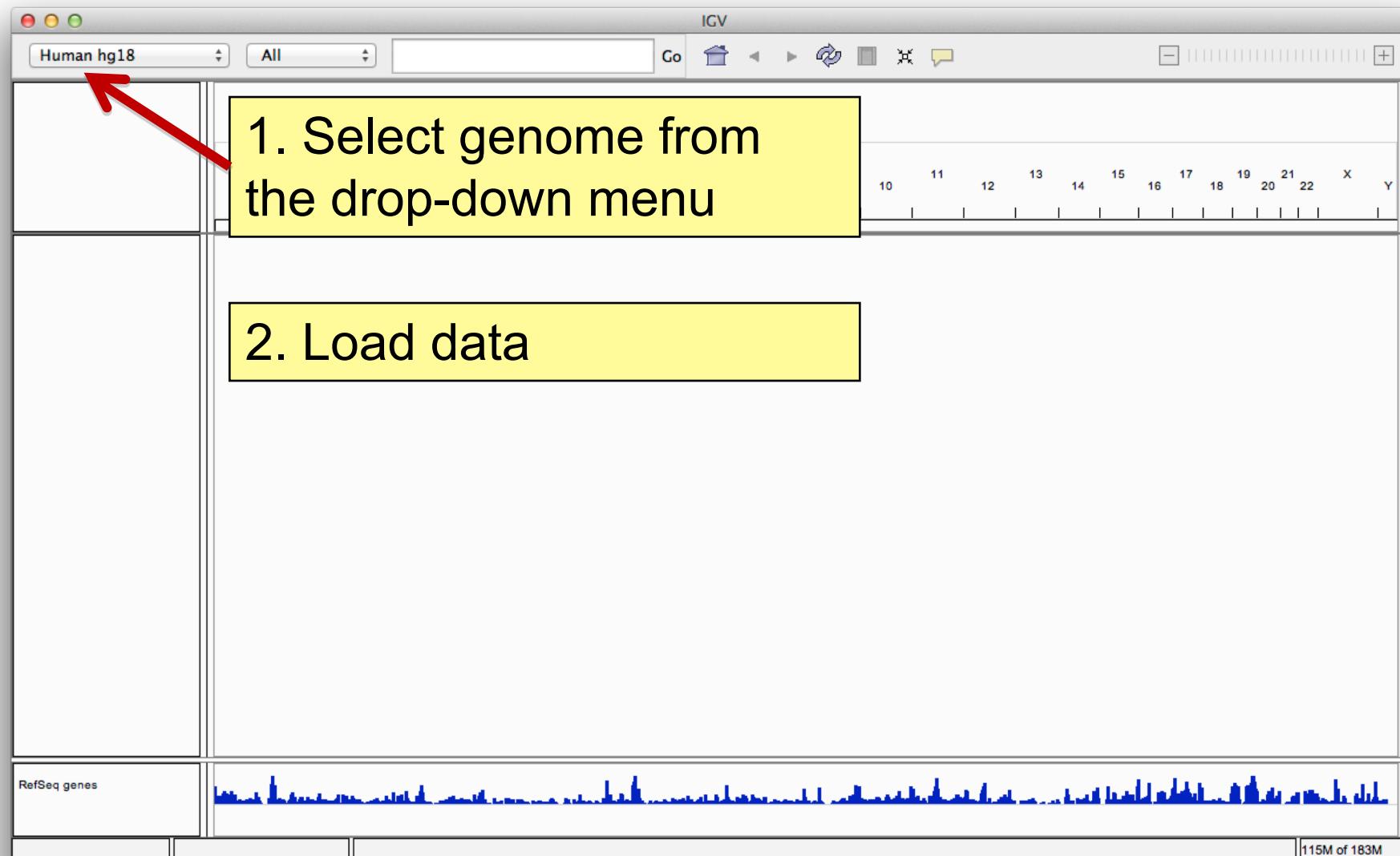


Load data

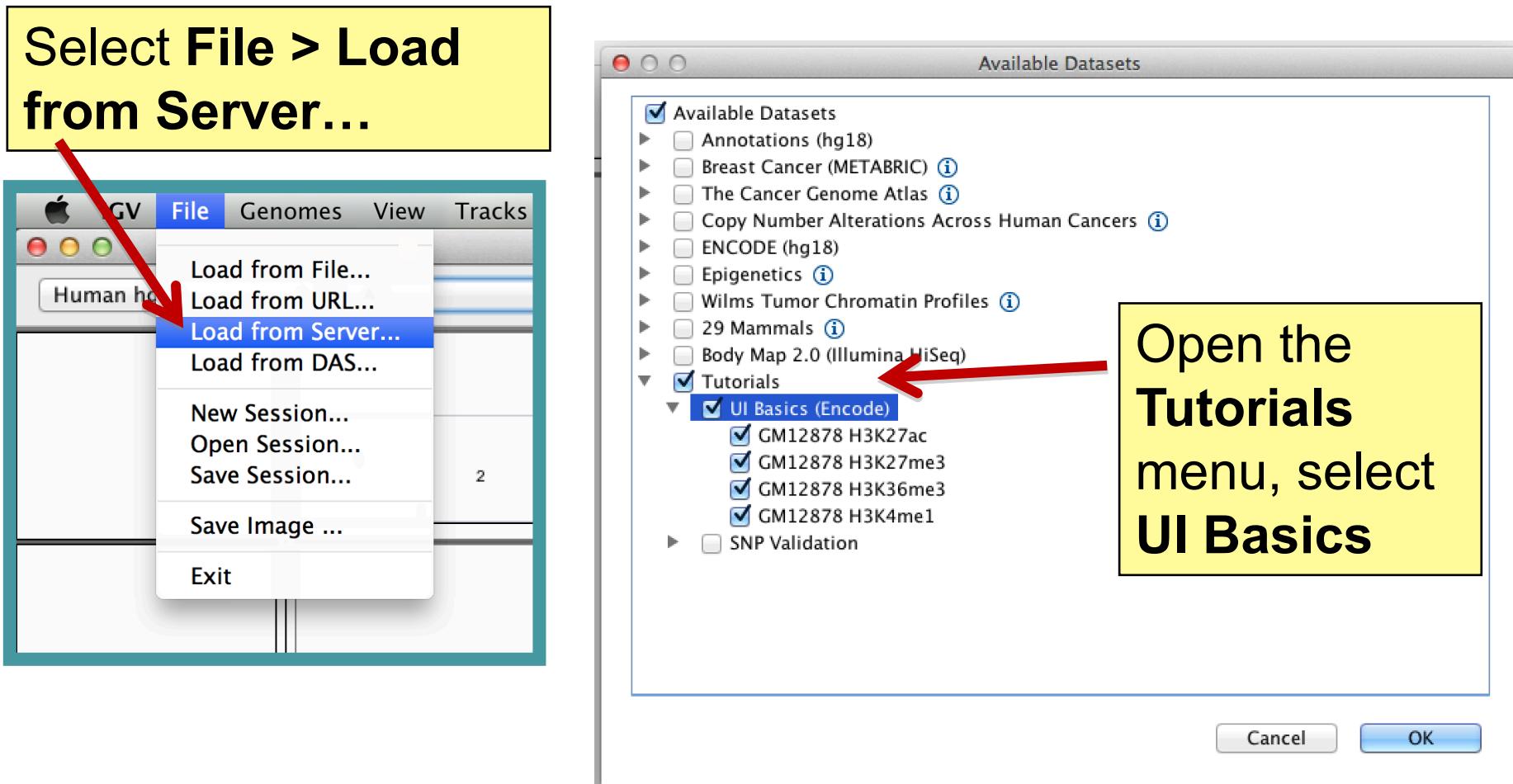

Screen layout

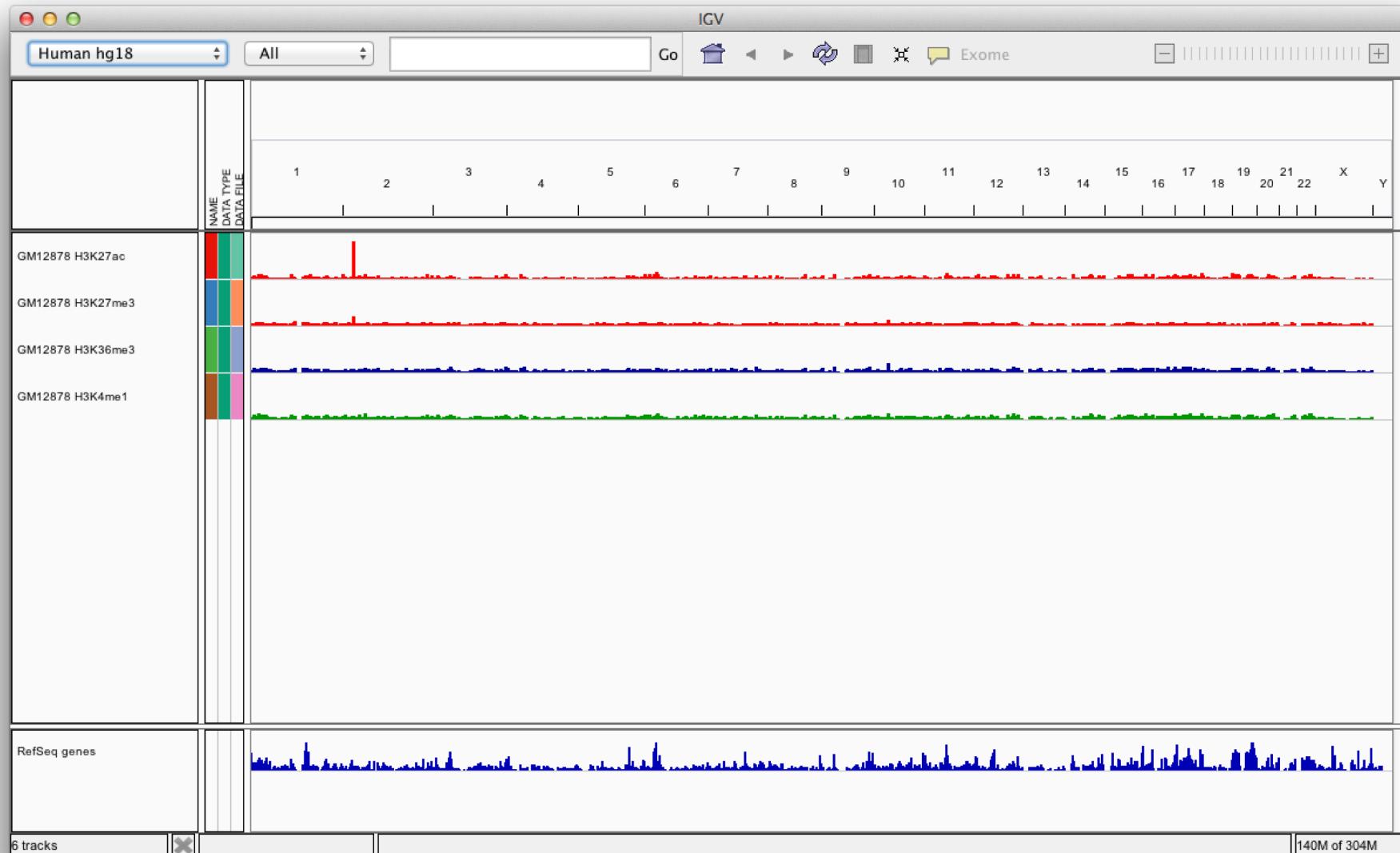

Screen layout

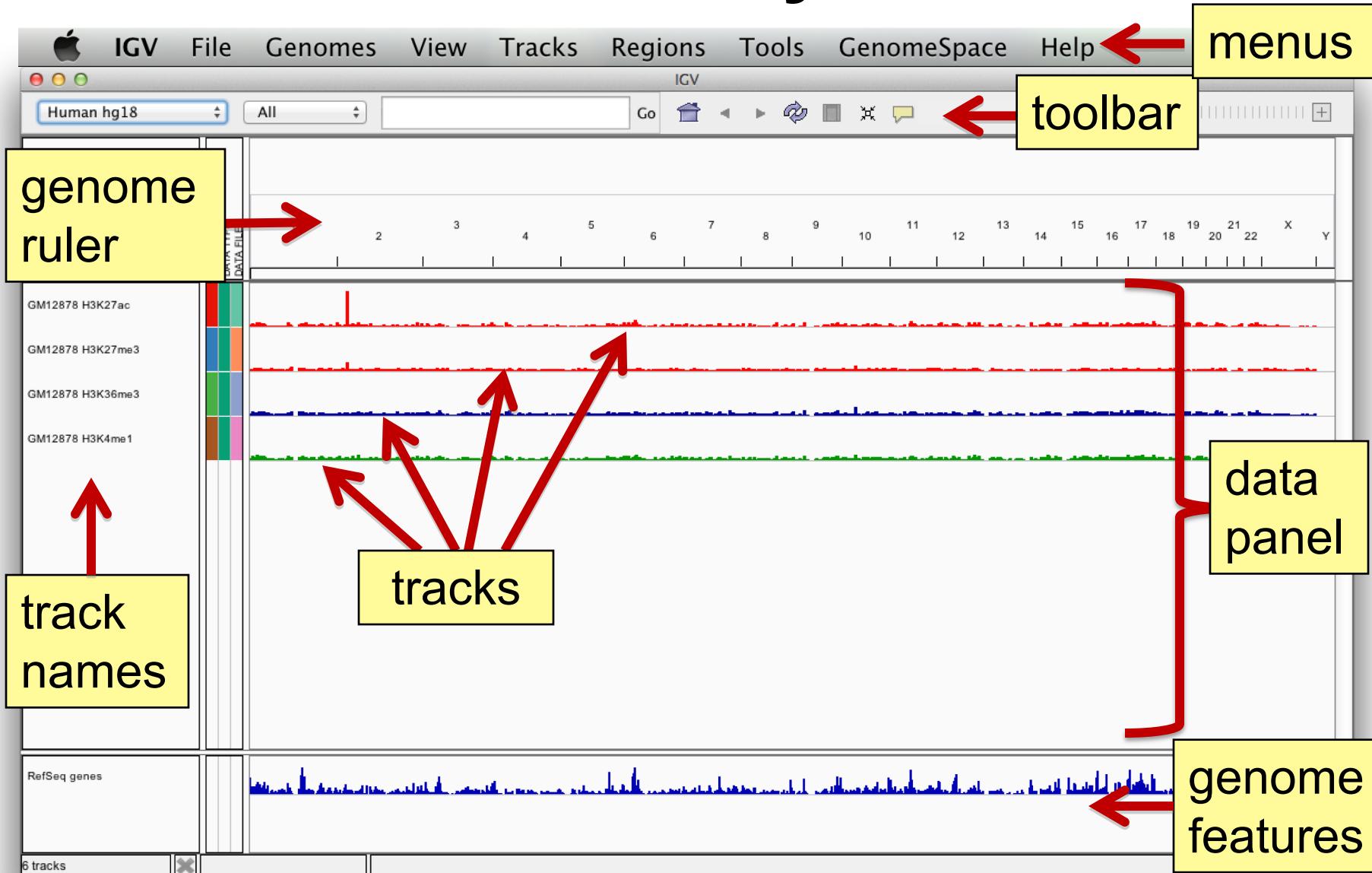

File formats and track types


- The **file format** defines the track type.
- The **track type** determines the display options

- [BAM](#)
- [BED](#)
- [BedGraph](#)
- [bigBed](#)
- [bigWig](#)
- [Birdsuite Files](#)
- [CBS](#)
- [CN](#)
- [Cufflinks Files](#)
- [Custom File Formats](#)
- [Cytoband](#)
- [FASTA](#)
- [GCT](#)
- [genePred](#)
- [GFF](#)
- [GISTIC](#)
- [Goby](#)
- [GWAS](#)
- [IGV](#)
- [LOH](#)
- [MAF](#)
- [Merged BAM File \(.bam.list\)](#)
- [MUT](#)
- [PSL](#)
- [RES](#)
- [SAM](#)
- [Sample Information](#)
- [SEG](#)
- [SNP](#)
- [TAB](#)
- [TDF](#)
- [Track Line](#)
- [Type Line](#)
- [VCF](#)
- [WIG](#)


- For current list see: www.broadinstitute.org/igv/FileFormats

Viewing alignments

Whole chromosome view

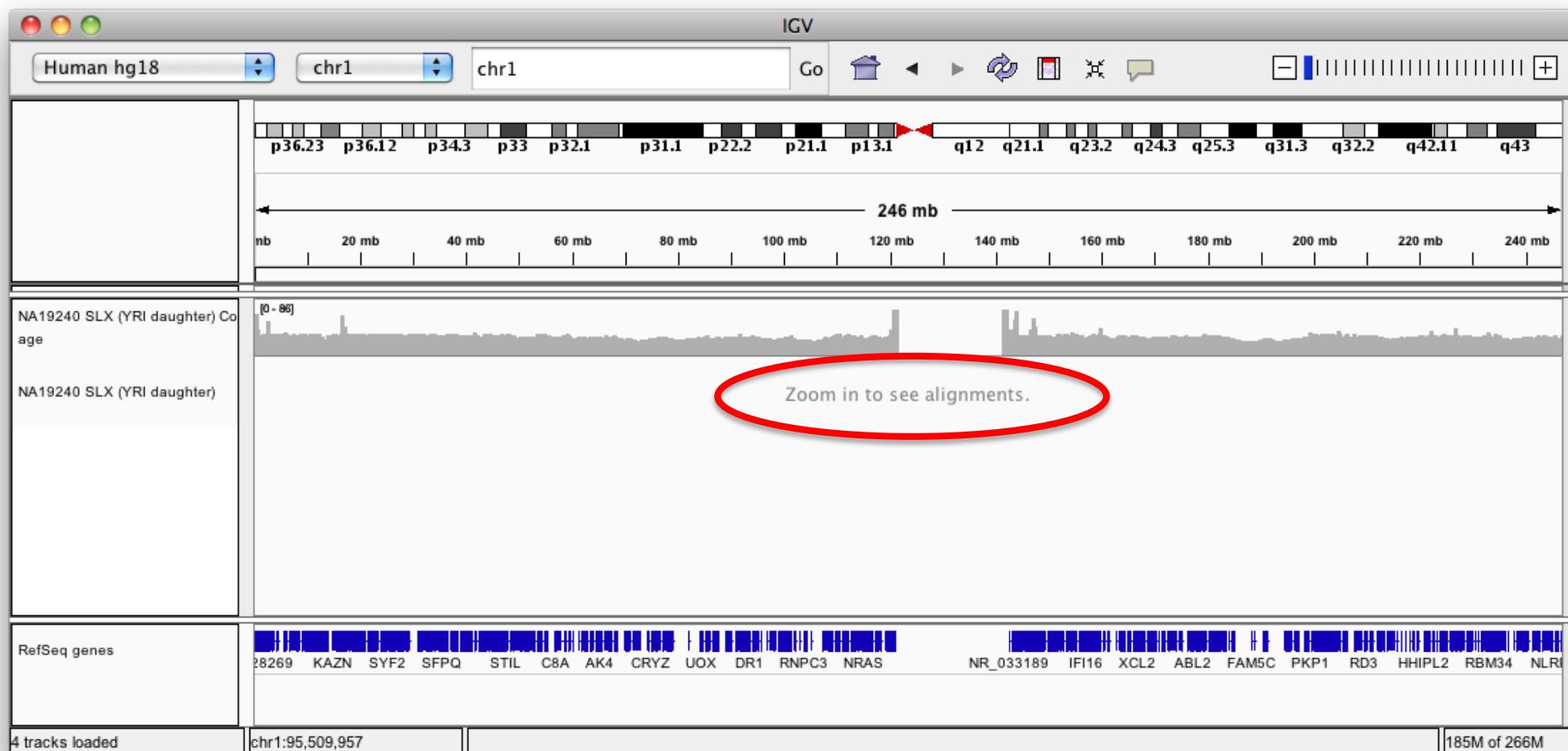

Viewing alignments – Zoom in

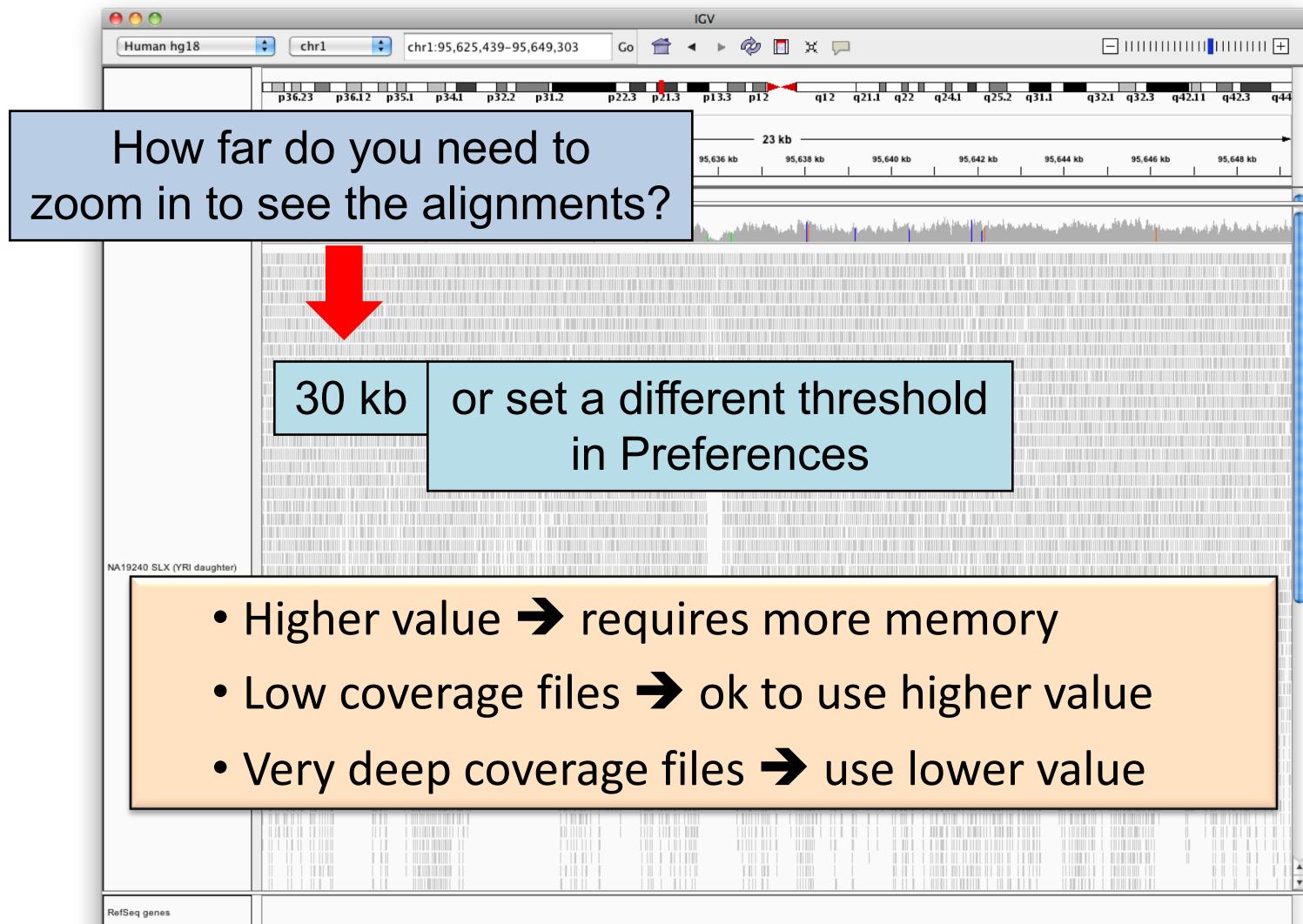


Viewing alignments – Zoom in



SNVs and Structural variations

- Important metrics for evaluating the validity of SNVs:
 - Coverage
 - Amount of support
 - Strand bias / PCR artifacts
 - Mapping qualities
 - Base qualities
- Important metrics for evaluating SVs:
 - Coverage
 - Insert size
 - Read pair orientation

Viewing SNPs and SNVs

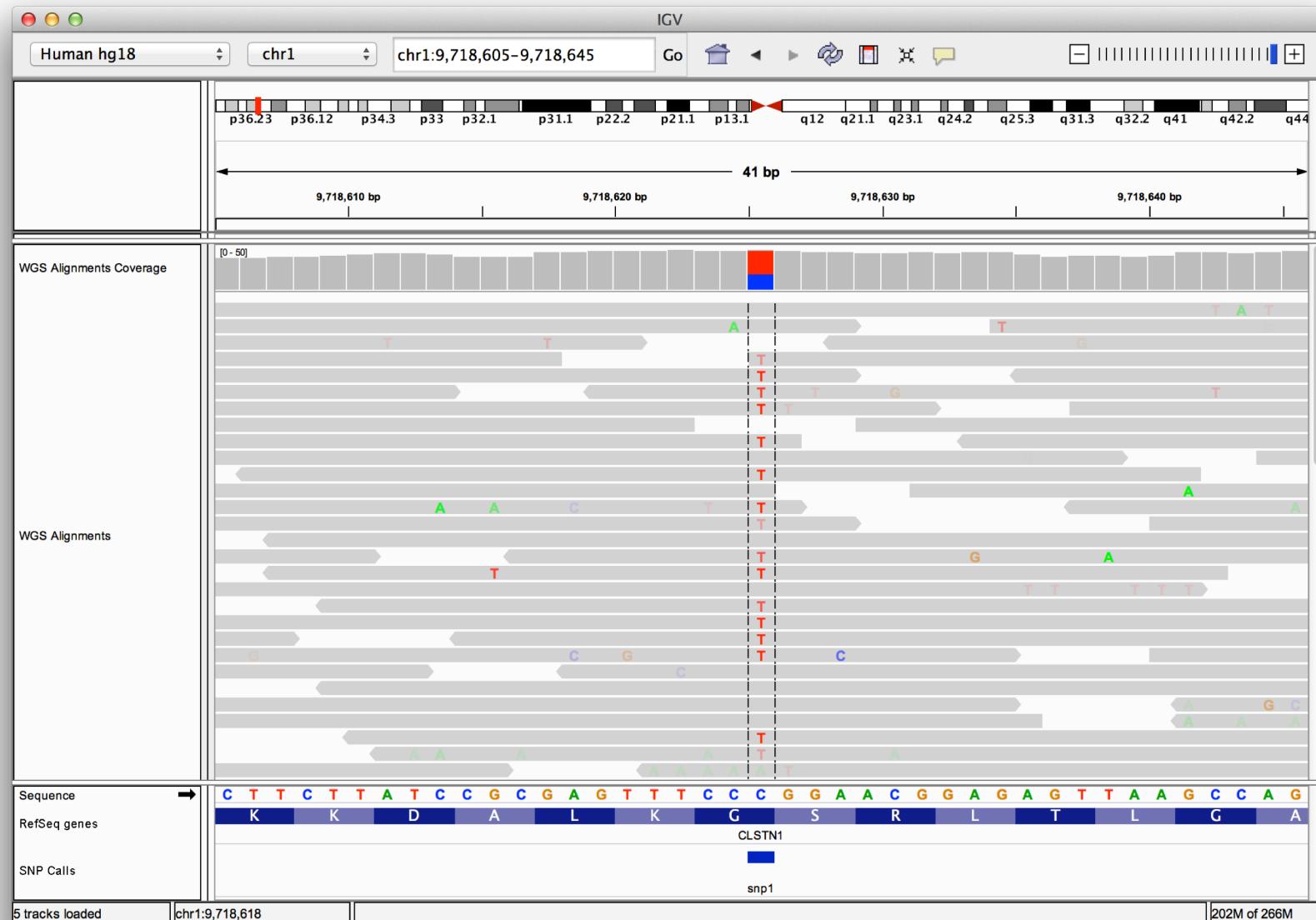


Viewing SNPs and SNVs




Viewing Structural Events

- Paired reads can yield evidence for genomic “structural events”, such as deletions, translocations, and inversions.
- Alignment coloring options help highlight these events based on:
 - Inferred insert size (template length)
 - Pair orientation (relative strand of pair)


Paired-end sequencing

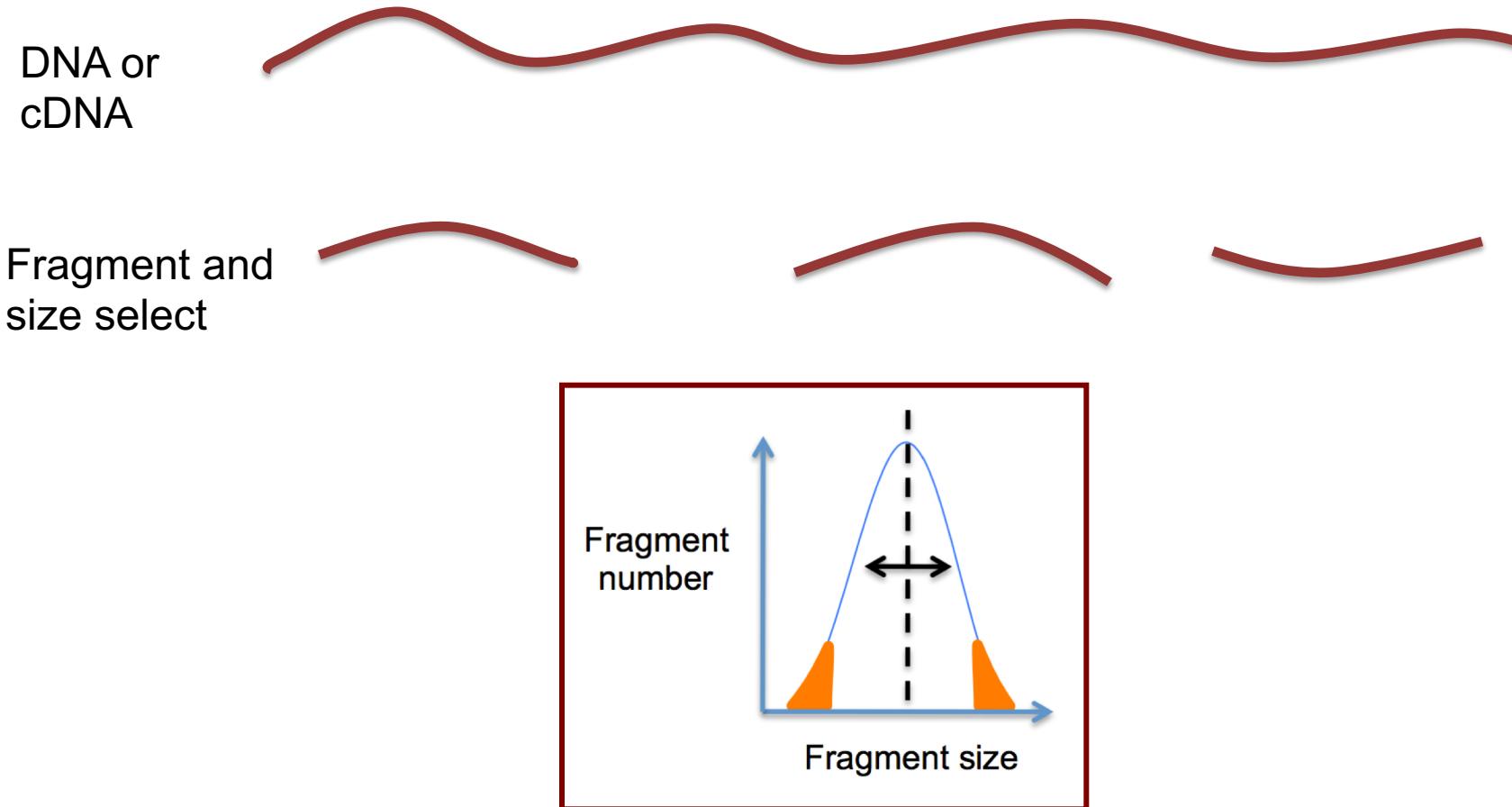

Paired-end sequencing

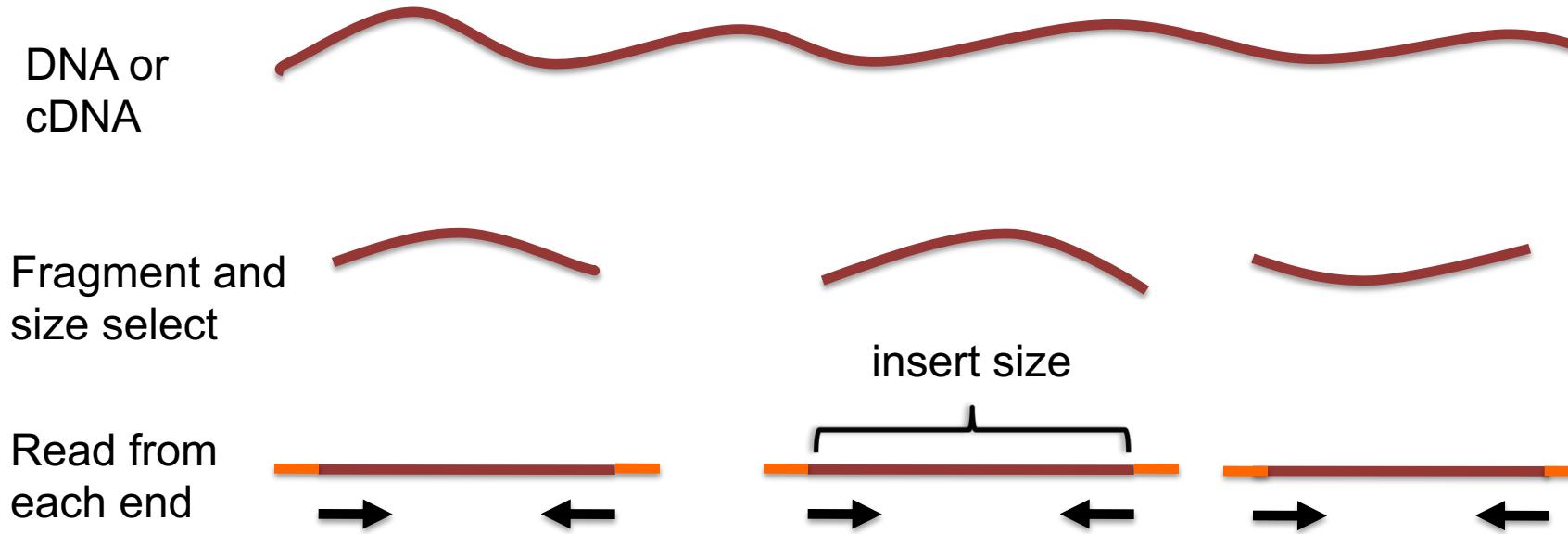

Paired-end sequencing

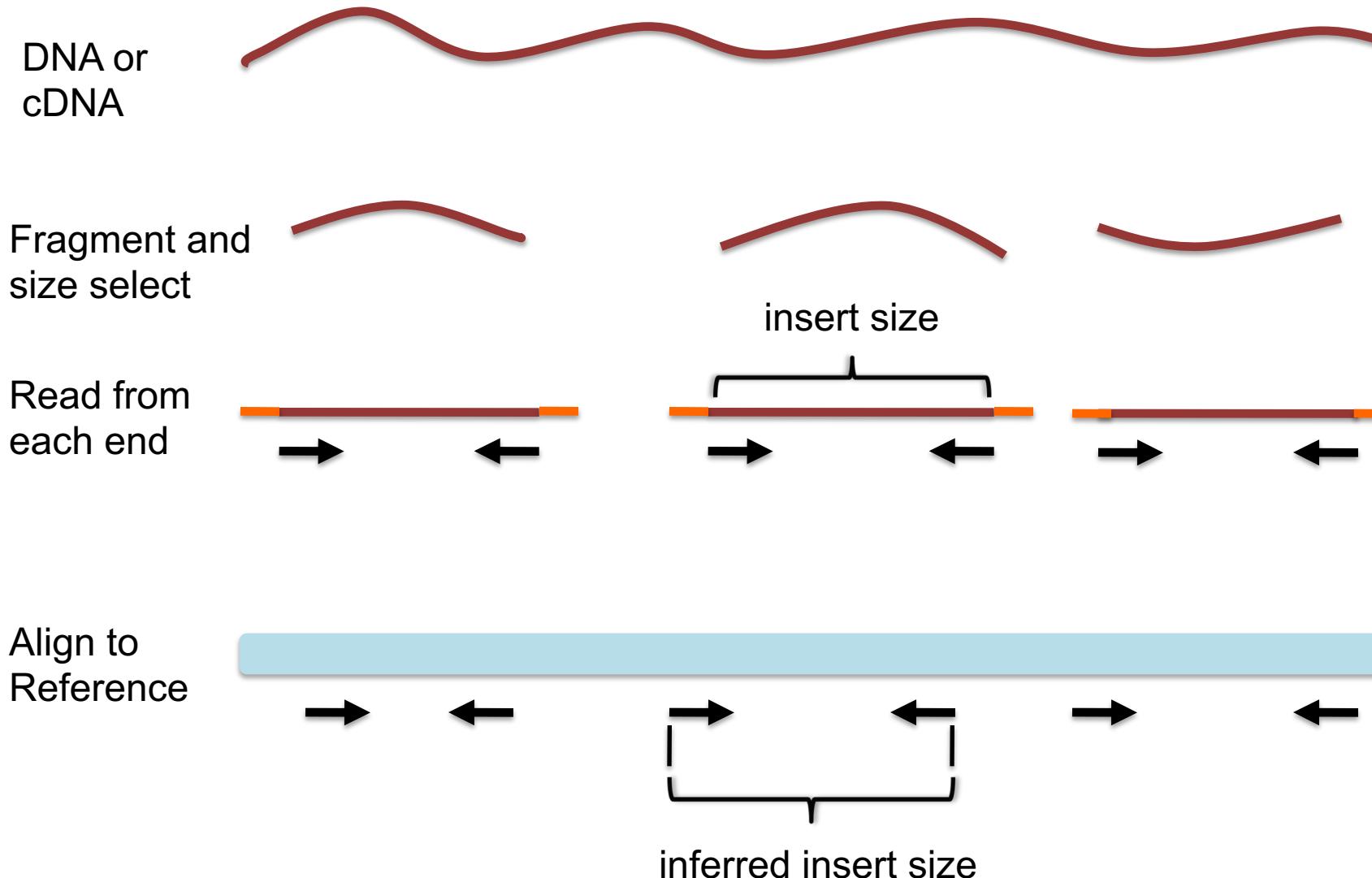

Interpreting inferred insert size

The “inferred insert size” can be used to detect structural variants including


- Deletions
- Insertions
- Inter-chromosomal rearrangements: (Undefined insert size)

Deletion


What is the effect of a deletion on inferred insert size?

Deletion

Reference
Genome

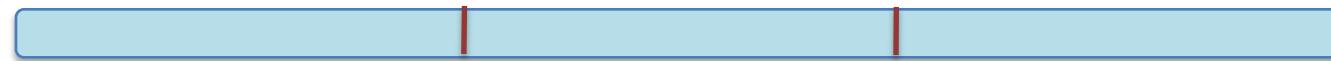


Subject

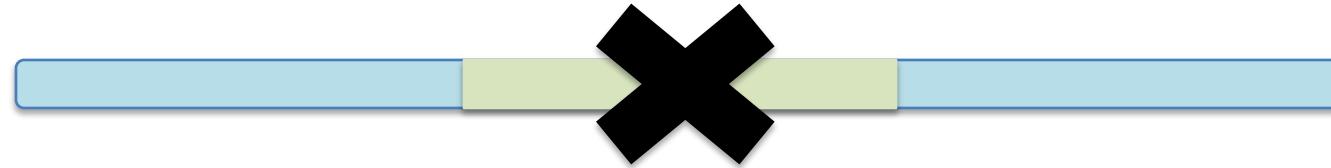

Deletion


Reference
Genome

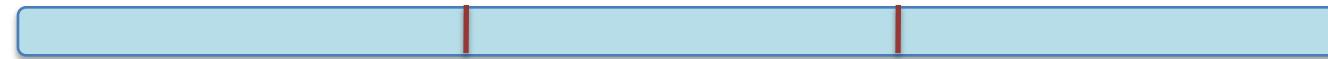

Subject




Deletion

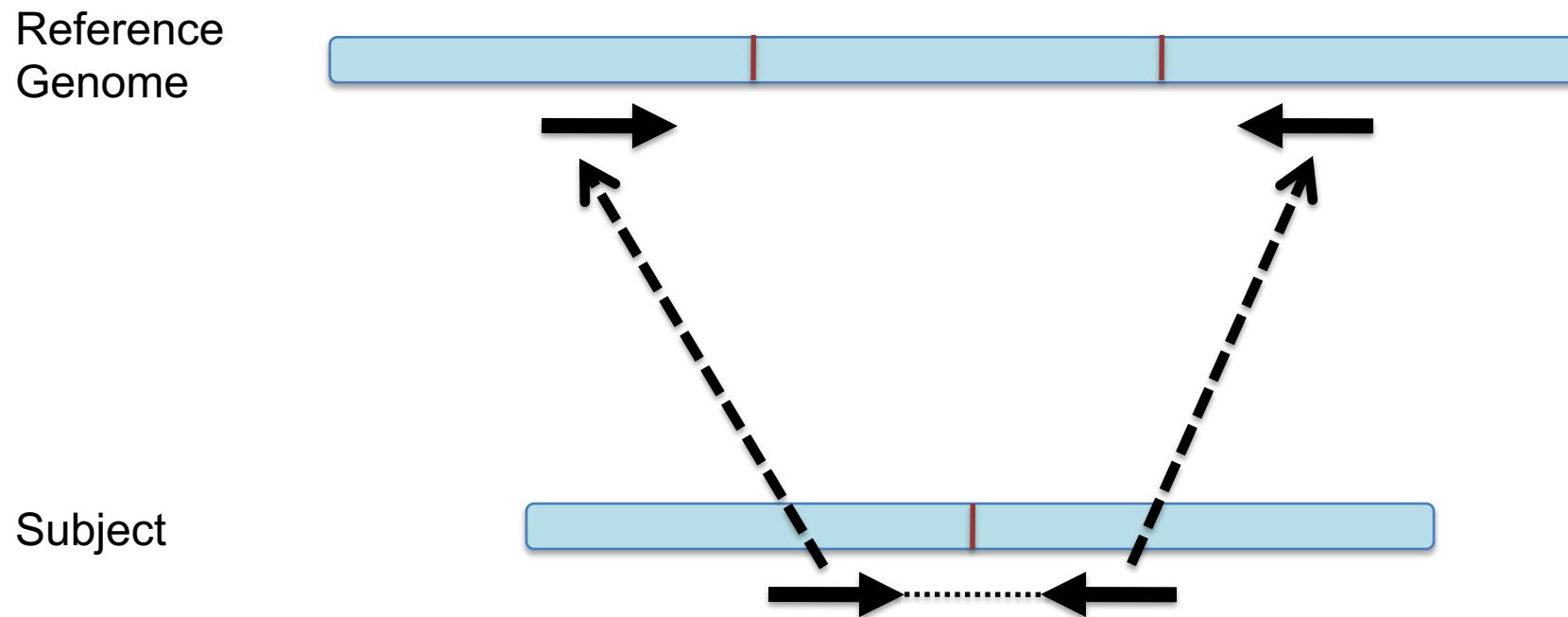

Deletion


Inferred insert size is > expected value

Reference
Genome

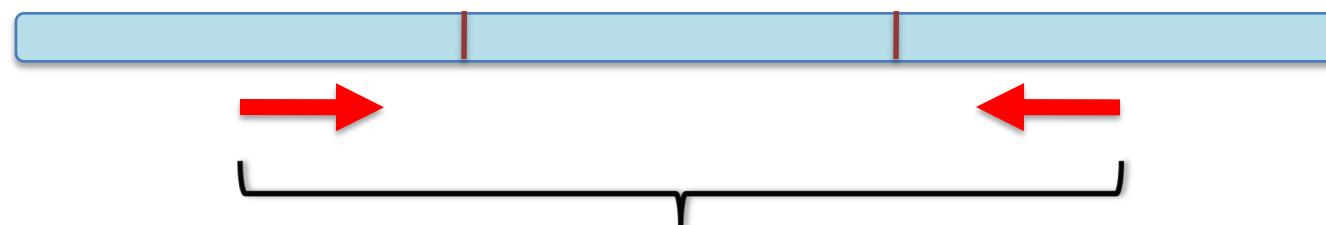


inferred insert size

Subject

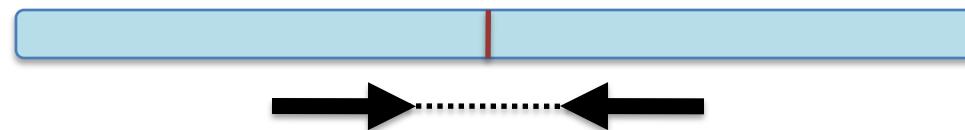


expected insert size

Color by insert size

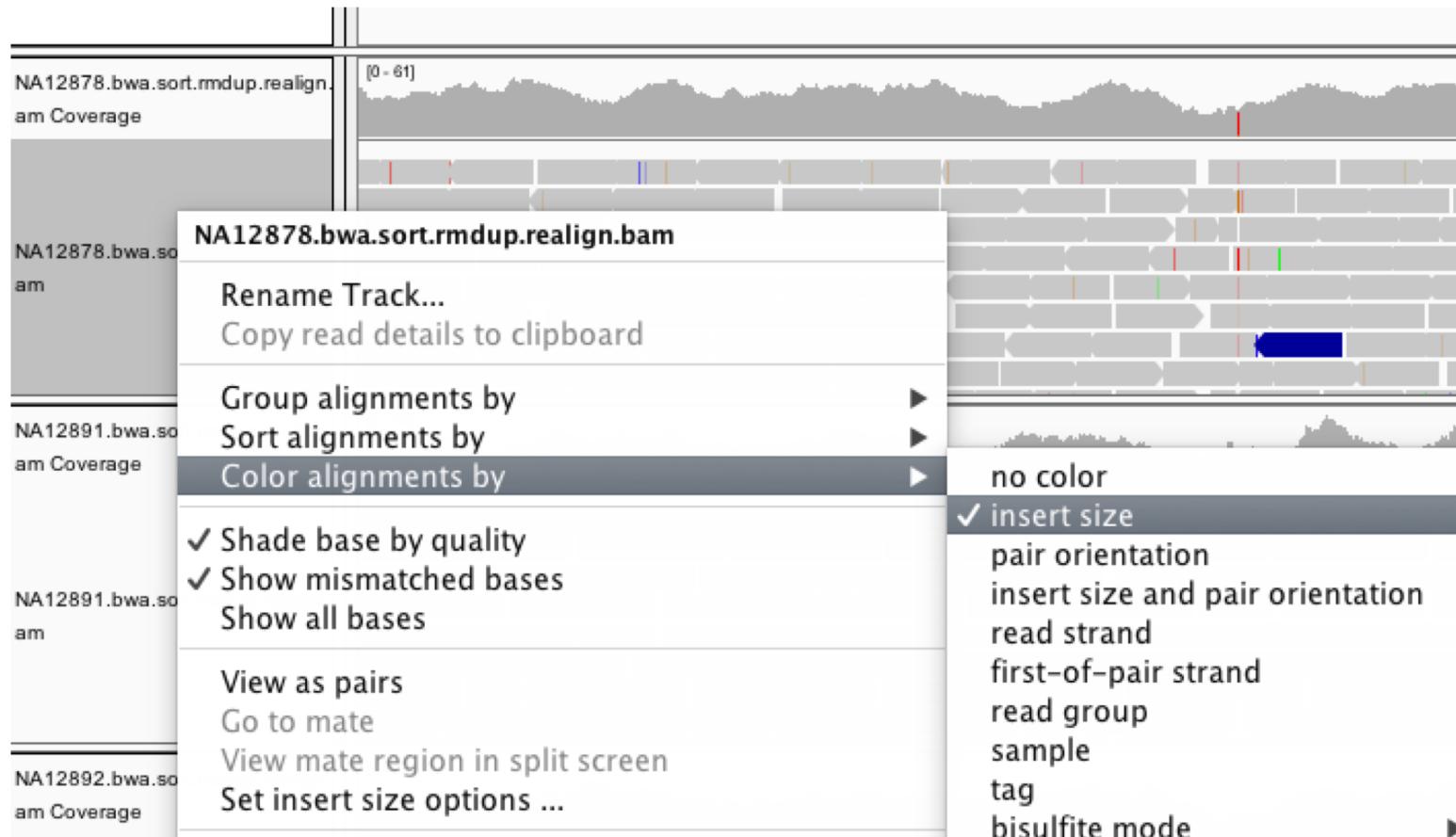


Deletion

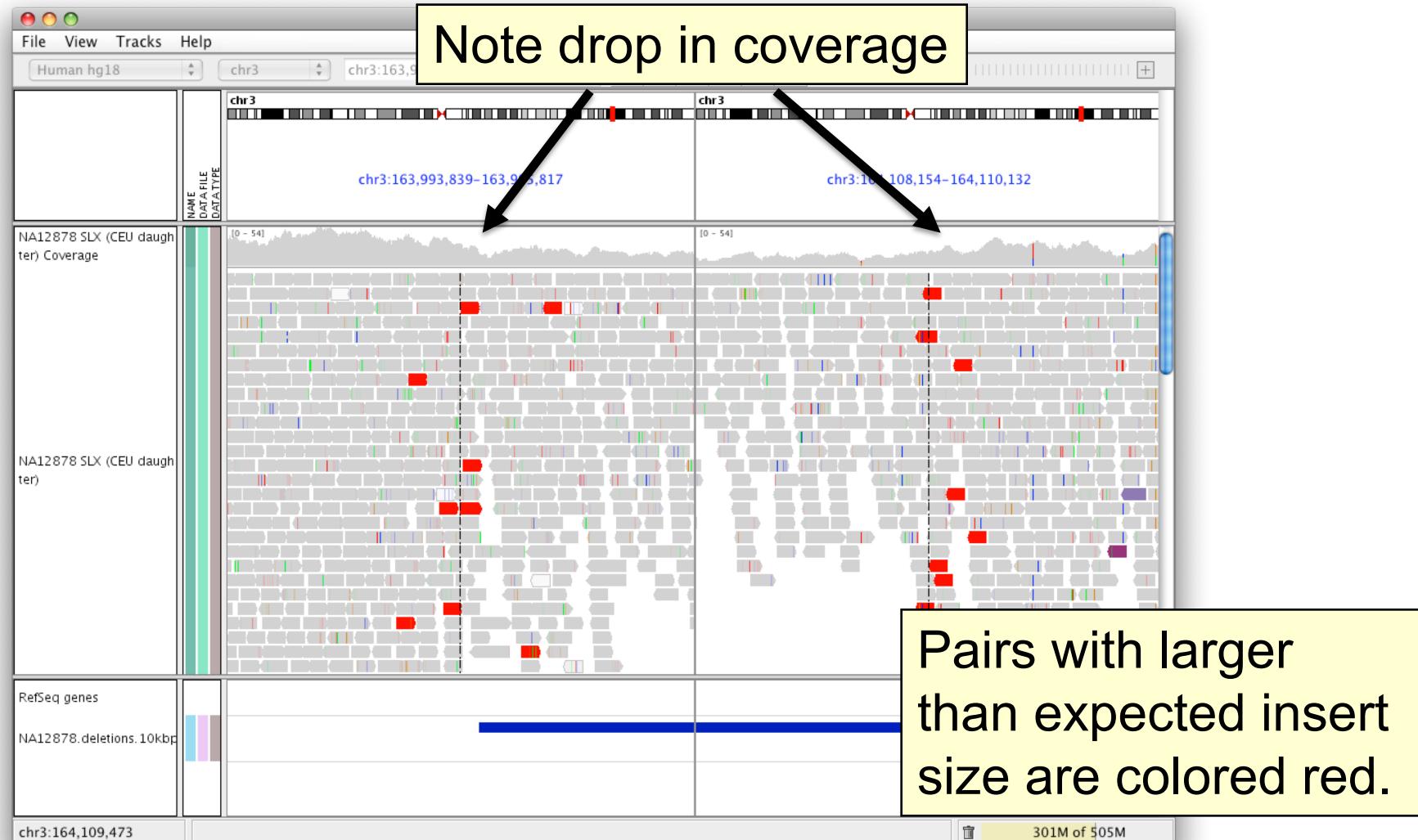


Insert size color scheme

- Smaller than expected insert size: 
- Larger than expected insert size: 
- Pairs on different chromosomes


Each end colored by chromosome of its mate

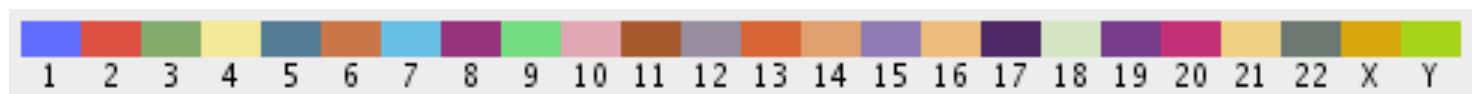

Rearrangement

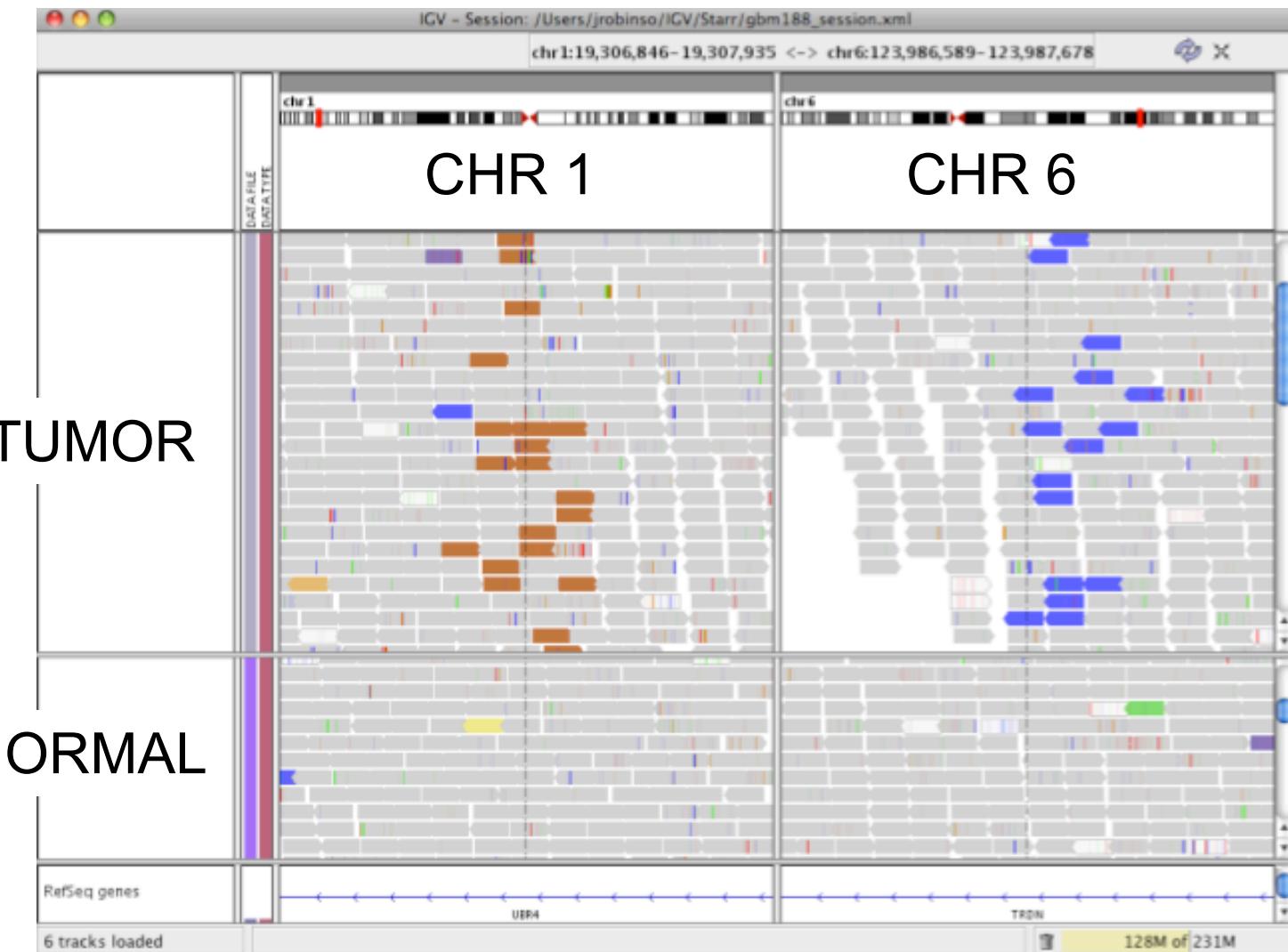

Rearrangement

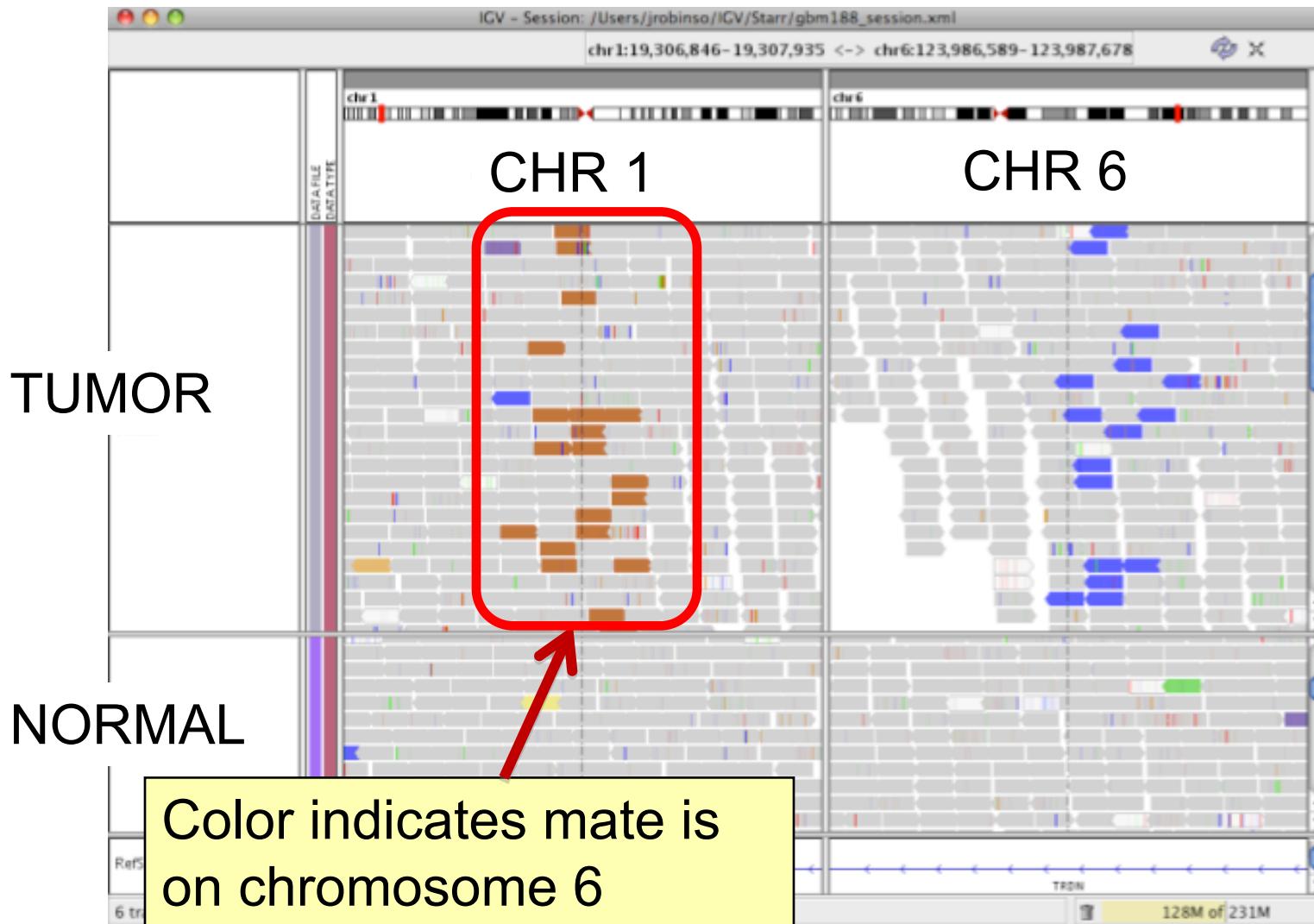

Interpreting Read-Pair Orientations

Orientation of paired reads can reveal structural events:

- Inversions
- Duplications
- Translocations
- Complex rearrangements

Orientation is defined in terms of

- read strand, left *vs* right, *and*
- read order, first *vs* second

Inversion


Reference
genome

Inversion

Reference
genome



Inversion

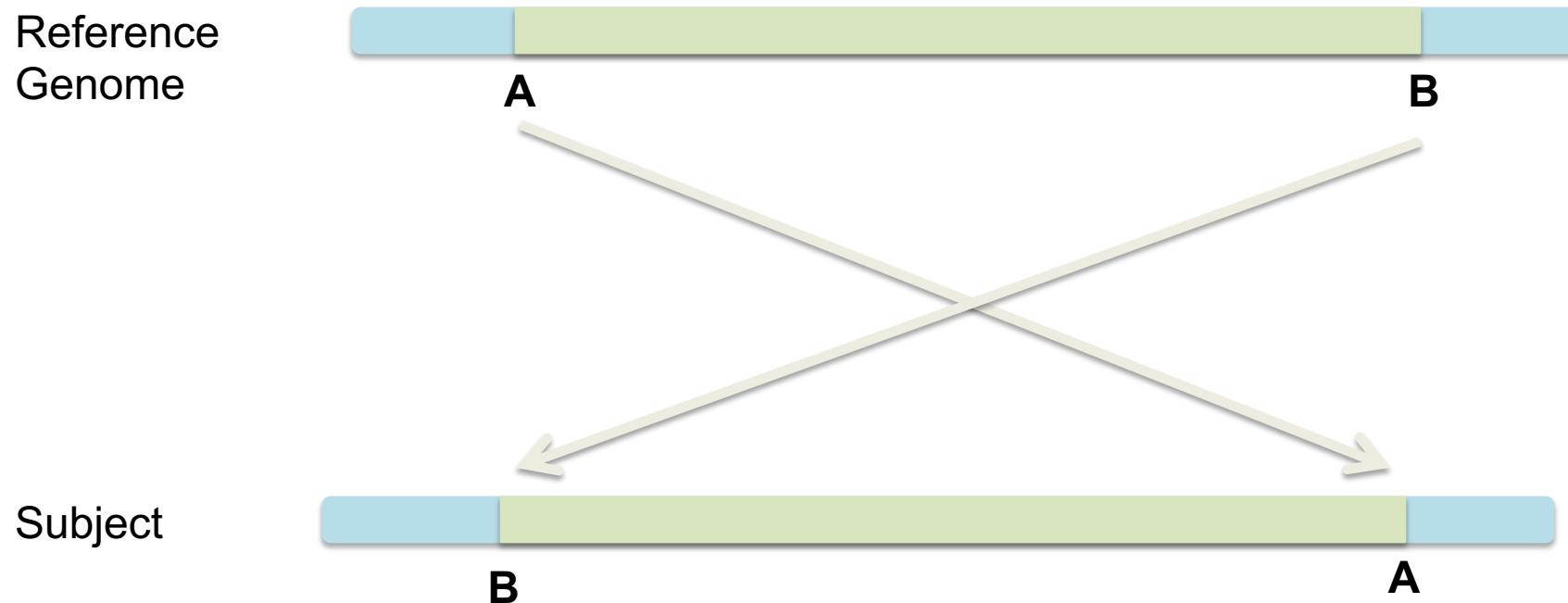

Inversion


Reference
Genome




Subject




Inversion

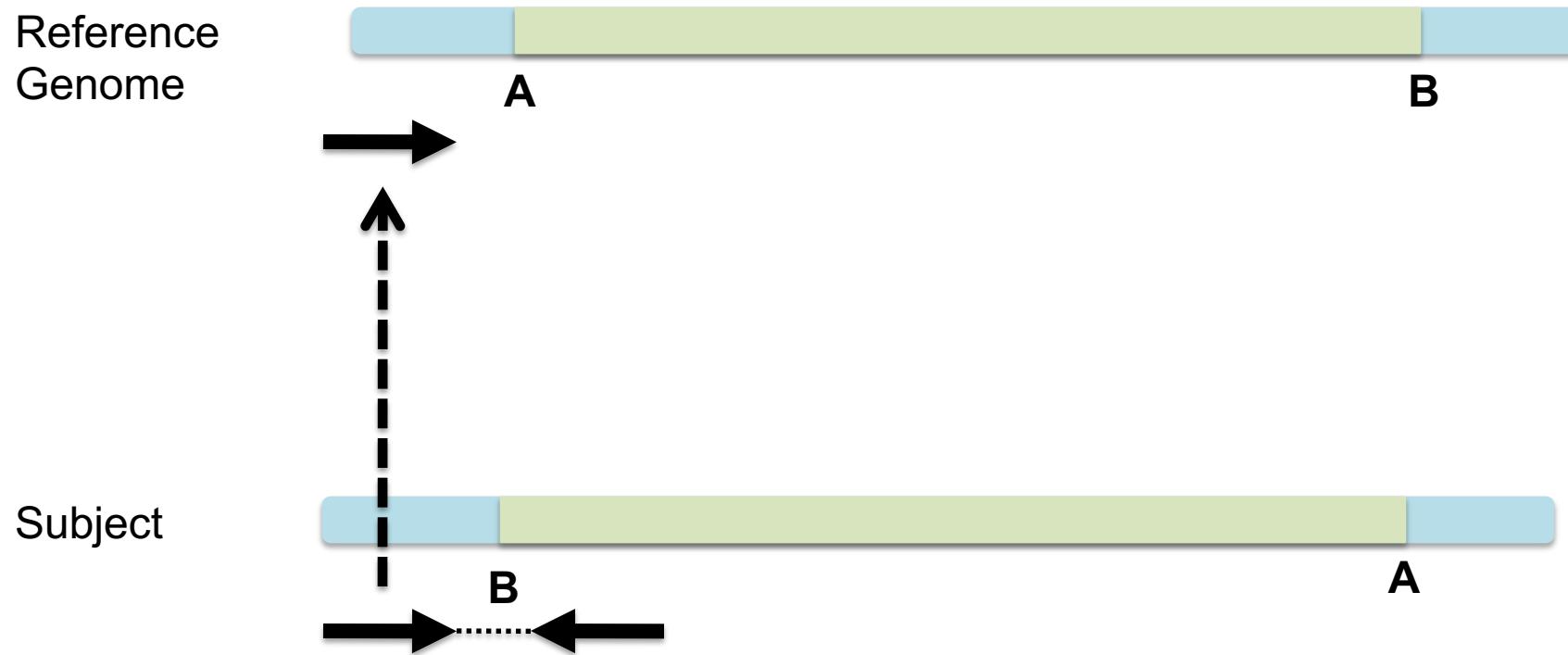


Inversion

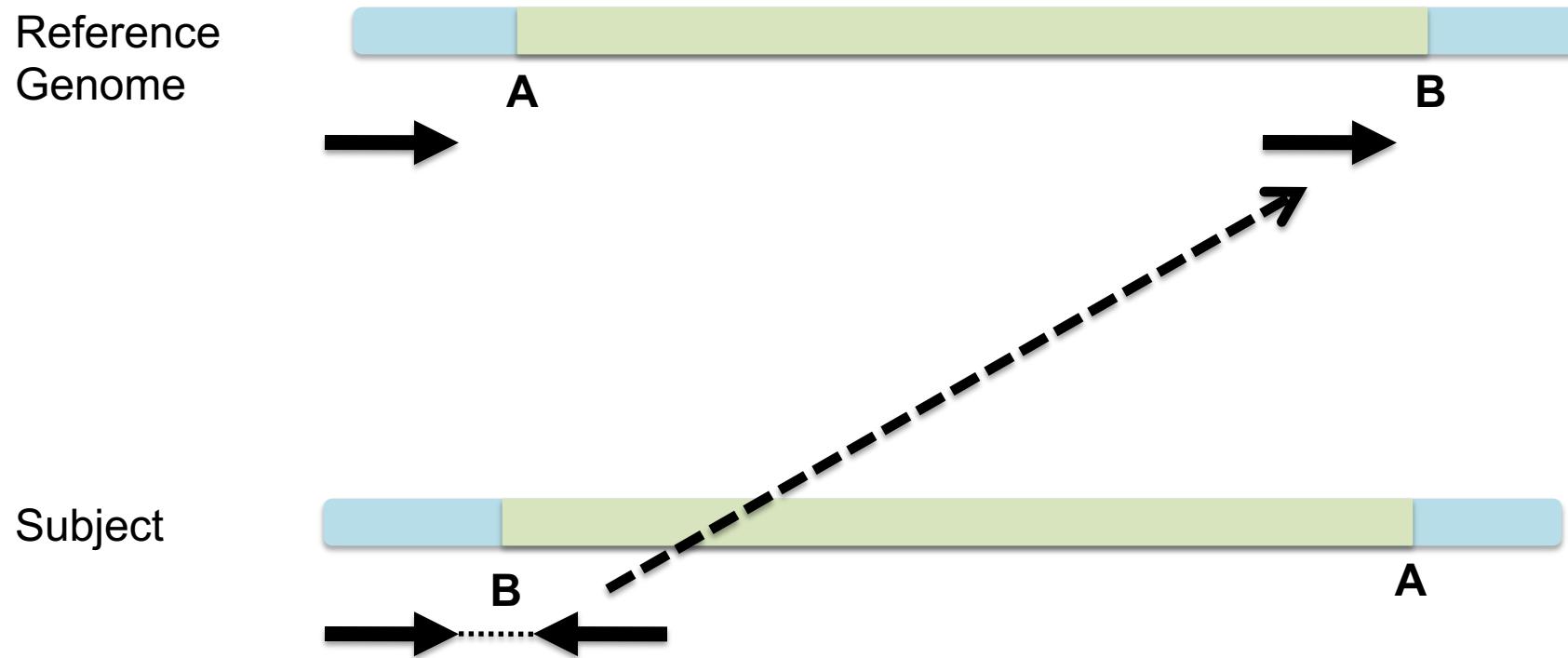


Inversion

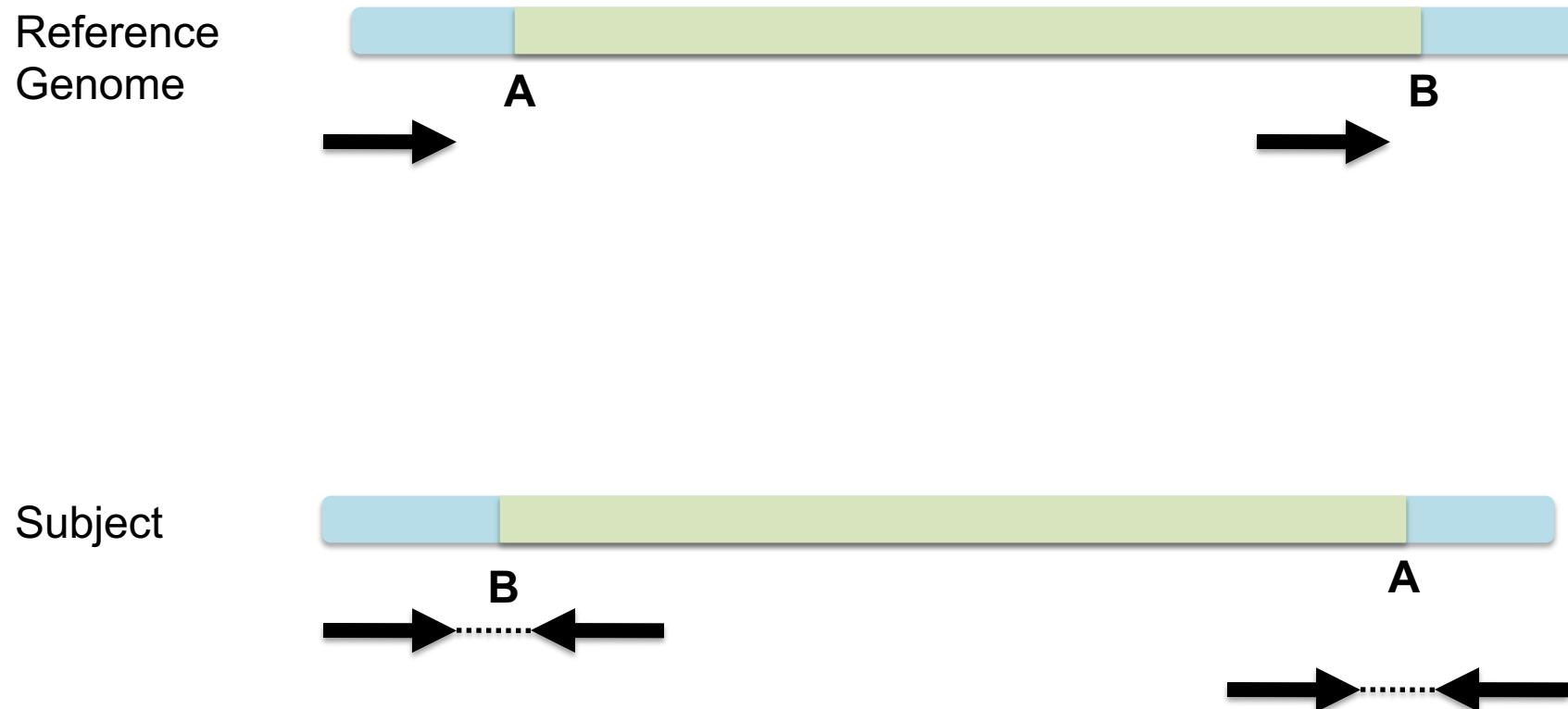


Inversion

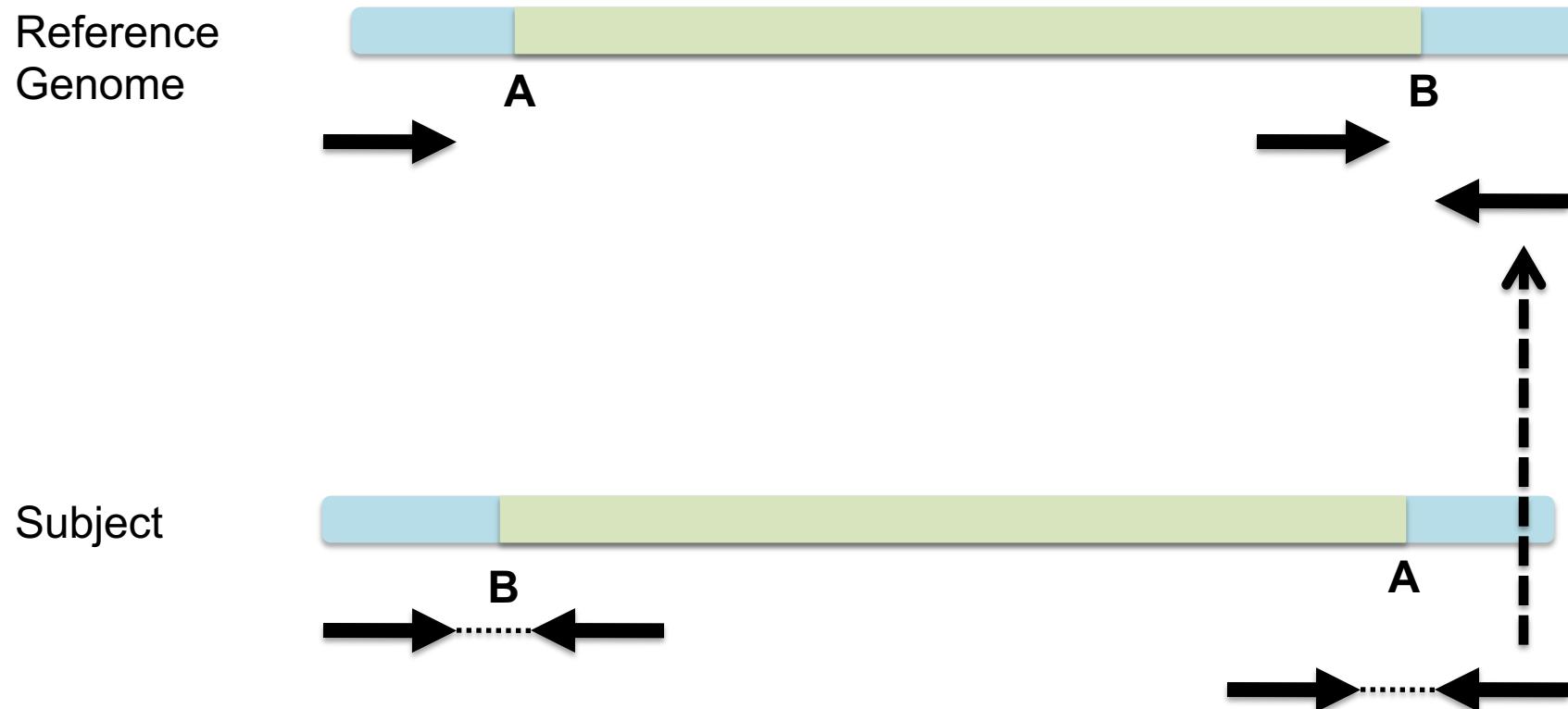

Inversion

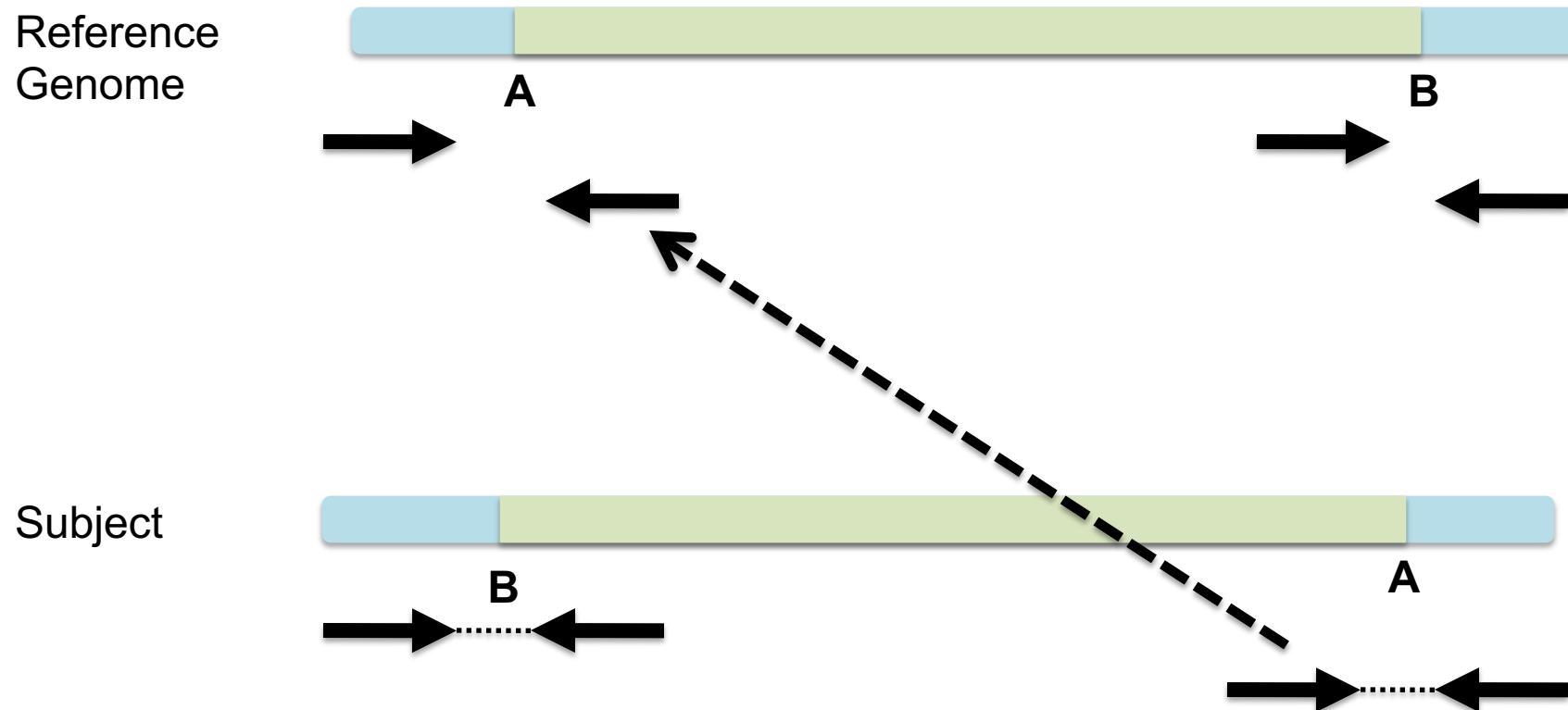

Inversion

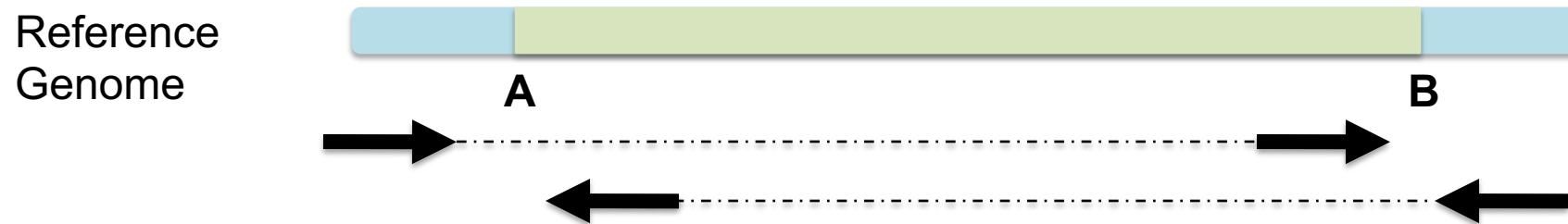


Inversion




Anomaly: expected orientation of pair is inward facing (→ ←)

Inversion

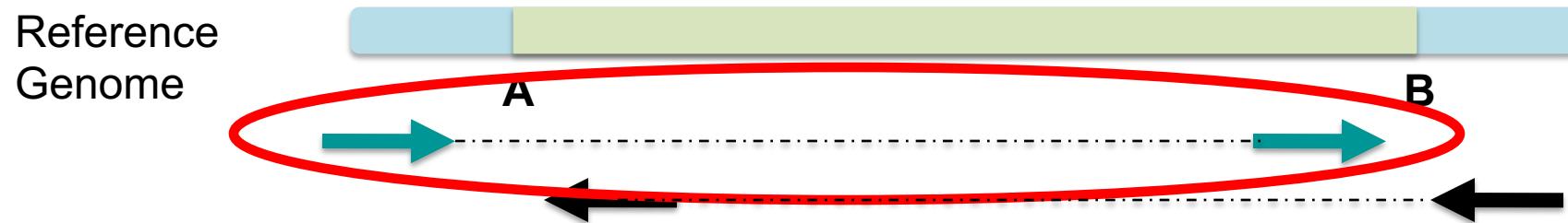


“Left” side pair


Inversion

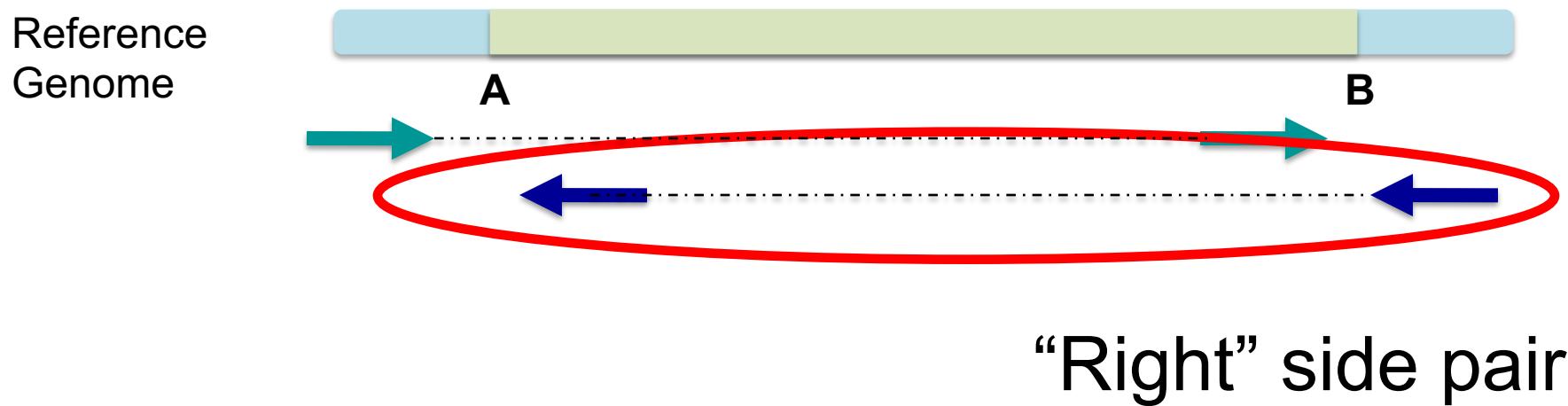


Color by pair orientation

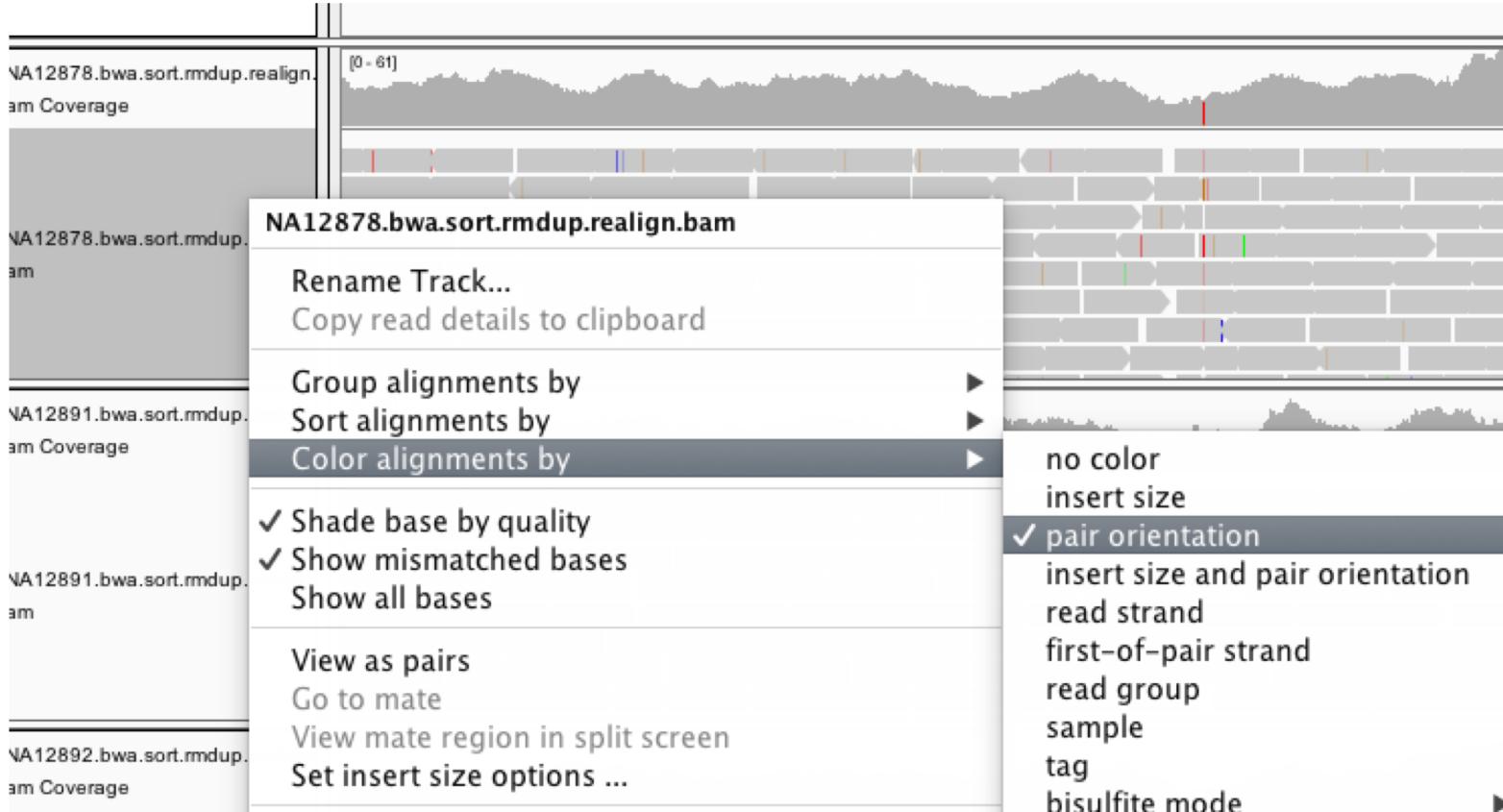

Inversion

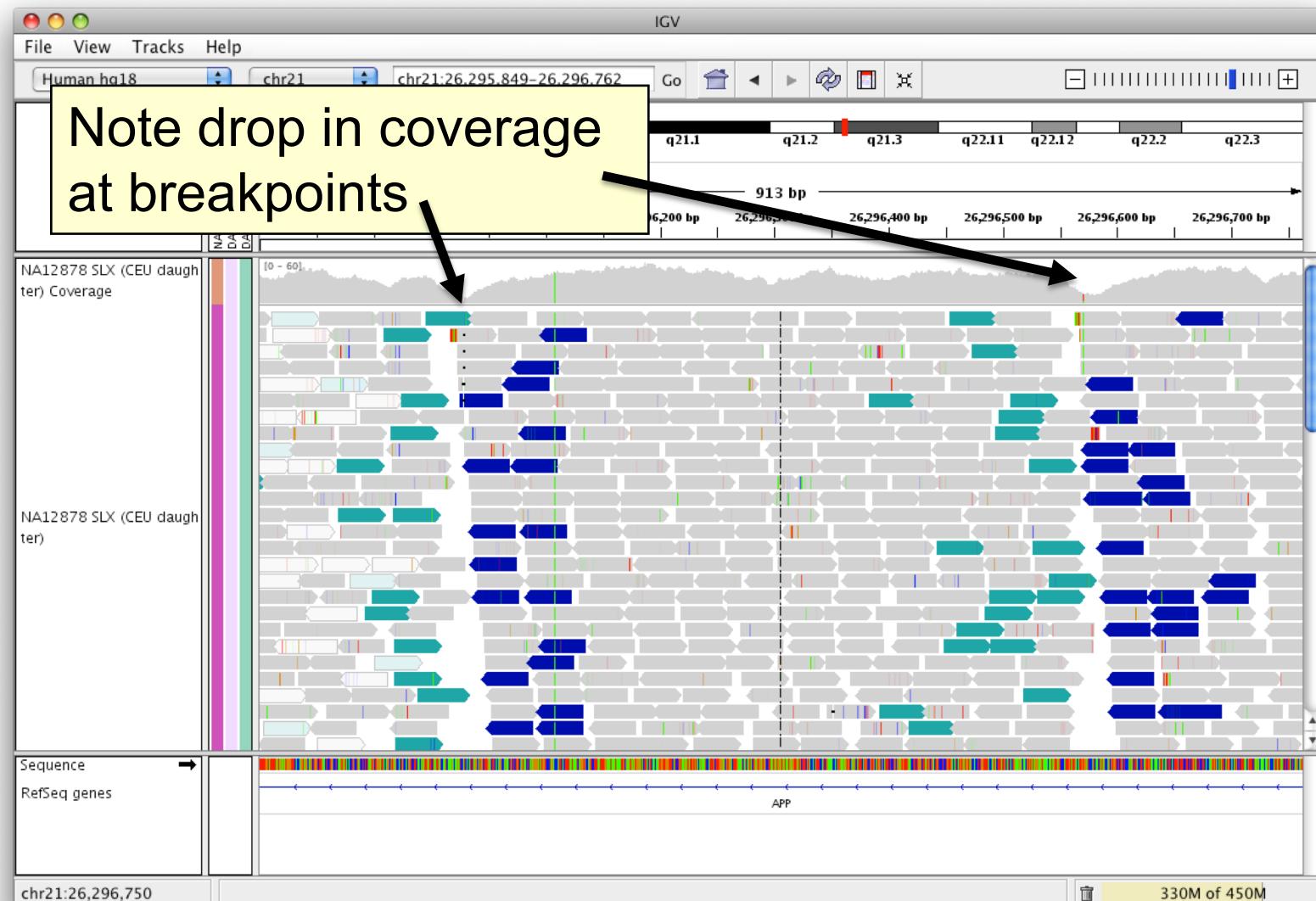

Interpretation of read pair orientations

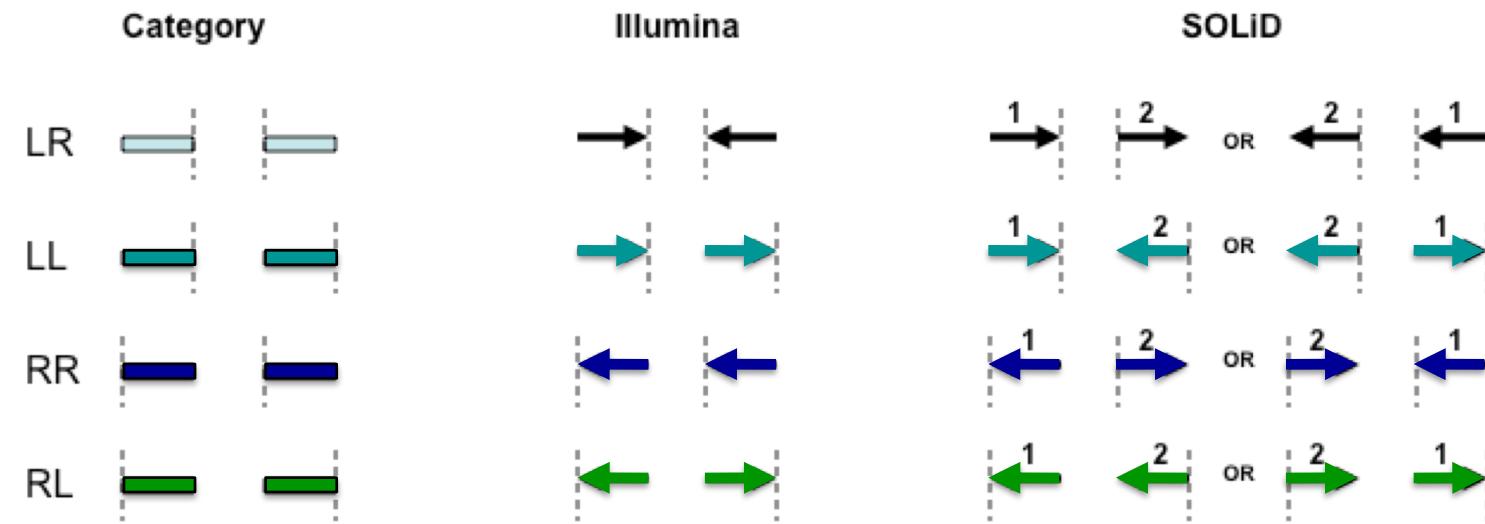

LR

Normal reads.

The reads are left and right (respectively) of the unsequenced part of the sequenced DNA fragment when aligned back to the reference genome.

LL,RR

Implies inversion in sequenced DNA with respect to reference.


RL

Implies duplication or translocation with respect to reference.

These categories only apply to reads where both mates map to the same chromosome.


Figure courtesy of Bob Handsaker

Long read considerations

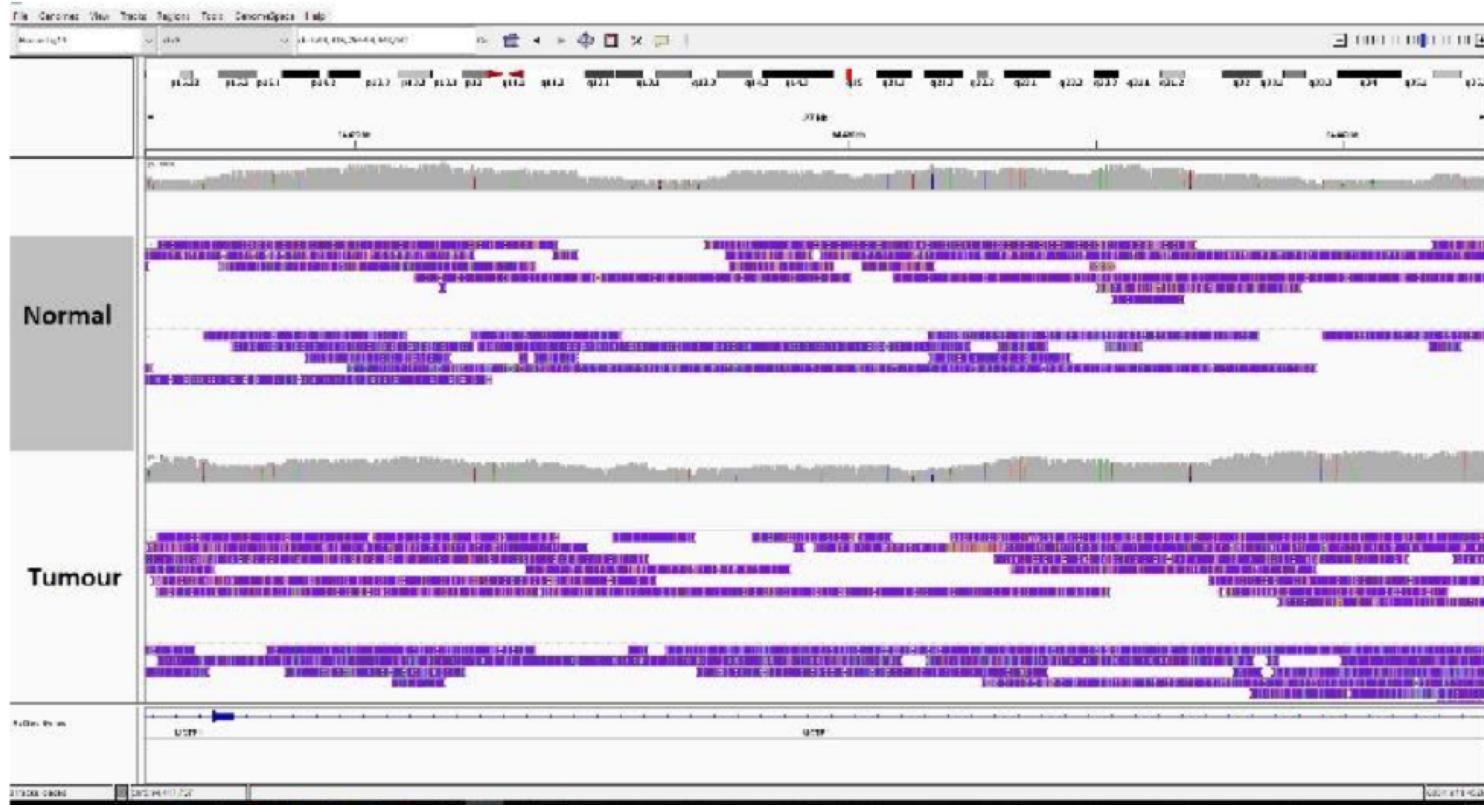

- Commonly see lots of small indels and single base errors that are simply noise
 - Can be removed to be able to view the data more cleanly

Long read considerations

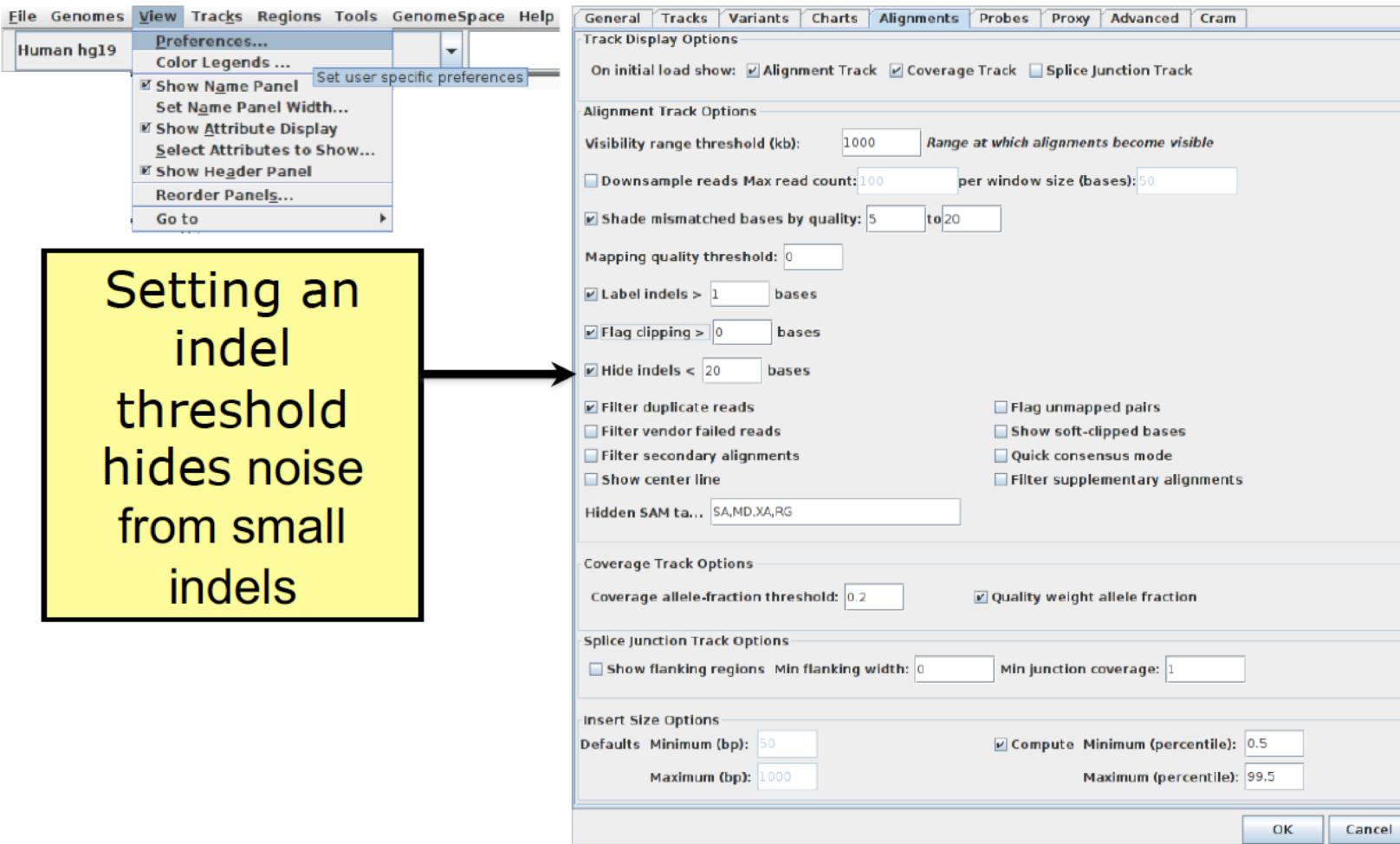

File Genomes View Tracks Regions Tools GenomeSpace Help

Human hg19

Human hg19

Preferences... Set user specific preferences

Show Name Panel
Show Attribute Display
Select Attributes to Show...
Show Header Panel
Reorder Panels...
Go to

Track Display Options

On initial load show: Alignment Track Coverage Track Splice Junction Track

Alignment Track Options

Visibility range threshold (kb): 1000 Range at which alignments become visible

Downsample reads Max read count: 100 per window size (bases): 50

Shade mismatched bases by quality: 5 to 20

Mapping quality threshold: 0

Label indels > 1 bases

Flag clipping > 0 bases

Hide indels < 20 bases

Filter duplicate reads

Filter vendor failed reads

Filter secondary alignments

Show center line

Flag unmapped pairs

Show soft-clipped bases

Quick consensus mode

Filter supplementary alignments

Hidden SAM ta... SA,MD,XA,RG

Coverage Track Options

Coverage allele-fraction threshold: 0.2 Quality weight allele fraction

Splice Junction Track Options

Show flanking regions Min flanking width: 0 Min junction coverage: 1

Insert Size Options


Defaults Minimum (bp): 50 Compute Minimum (percentile): 0.5

Maximum (bp): 1000 Maximum (percentile): 99.5

OK Cancel

Setting an indel threshold hides noise from small indels

Long read considerations

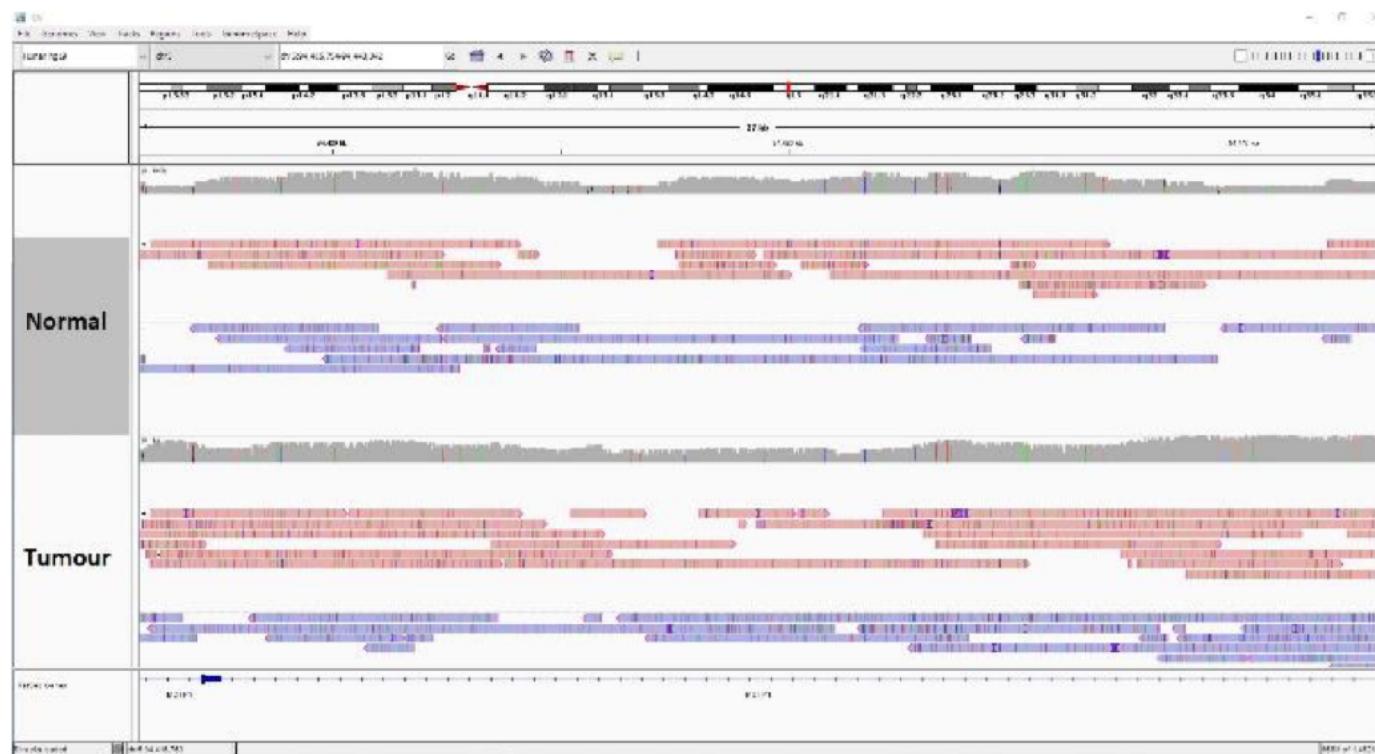

- Reads are not all purple dashes
- Next step would be to call a consensus at each position

Long read considerations

The screenshot shows the IGV software interface with the 'Track Display Options' dialog open. The 'Alignment Track Options' section is visible, containing various settings for alignment tracks. A yellow callout box highlights the 'Quick consensus mode' checkbox, which is checked. Other options shown include 'Filter duplicate reads', 'Filter vendor failed reads', 'Filter secondary alignments', 'Show center line', 'Flag unmapped pairs', 'Show soft-clipped bases', 'Filter supplementary alignments', and 'Hidden SAM ta... SA,MD,XA,RG'. The 'OK' and 'Cancel' buttons are at the bottom right of the dialog.

File Genomes View Tracks Regions Tools GenomeSpace Help

Human hg19

Preferences... Color Legends ... Set user specific preferences

Show Name Panel Set Name Panel Width... Show Attribute Display Select Attributes to Show... Show Header Panel Reorder Panels... Go to

General Tracks Variants Charts Alignments Probes Proxy Advanced Cram

Track Display Options

On initial load show: Alignment Track Coverage Track Splice Junction Track

Alignment Track Options

Visibility range threshold (kb): 1000 Range at which alignments become visible

Downsample reads Max read count: 100 per window size (bases): 50

Shade mismatched bases by quality: 5 to 20

Mapping quality threshold: 0

Label indels > 1 bases

Flag clipping > 0 bases

Hide indels < 20 bases

Filter duplicate reads

Filter vendor failed reads

Filter secondary alignments

Show center line

Flag unmapped pairs

Show soft-clipped bases

Quick consensus mode

Filter supplementary alignments

Hidden SAM ta... SA,MD,XA,RG

Coverage Track Options

Coverage allele-fraction threshold: 0.2 Quality weight allele fraction

Splice Junction Track Options

Show flanking regions Min flanking width: 0 Min junction coverage: 1

Insert Size Options


Defaults Minimum (bp): 50 Compute Minimum (percentile): 0.5

Maximum (bp): 1000 Maximum (percentile): 99.5

OK Cancel

Option for generating consensus sequences

Long read considerations

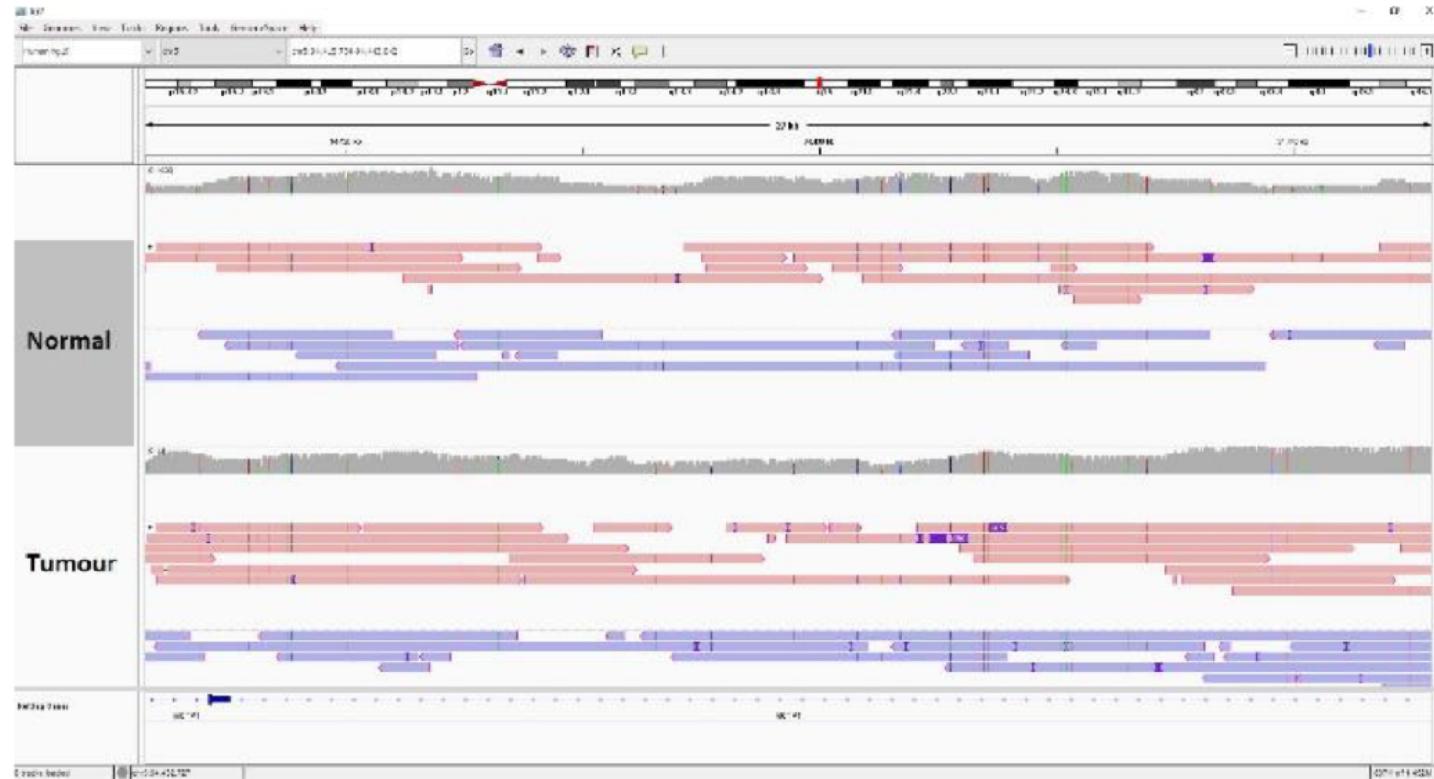

- Much easier to parse through the genomic data
- Large insertions and deletions are also labelled now

Manual Review Standard Operating Procedure (SOP) paper

© American College of Medical Genetics and Genomics

ARTICLE

Genetics
inMedicine

Open

Standard operating procedure for somatic variant refinement of sequencing data with paired tumor and normal samples

Erica K. Barnell, BS¹, Peter Ronning, BS¹, Katie M. Campbell, BS¹, Kilannin Krysiak, PhD^{1,2}, Benjamin J. Ainscough, PhD^{1,3}, Lana M. Sheta¹, Shahil P. Pema¹, Alina D. Schmidt, BS¹, Megan Richters, BS¹, Kelsy C. Cotto, BS¹, Arpad M. Danos, PhD¹, Cody Ramirez, BS¹, Zachary L. Skidmore, MEng¹, Nicholas C. Spies, BS¹, Jasreet Hundal, MS¹, Malik S. Sediqzad¹, Jason Kunisaki, BS¹, Felicia Gomez, PhD¹, Lee Trani, BS¹, Matthew Matlock, BS¹, Alex H. Wagner, PhD¹, S. Joshua Swamidass, MD/PhD^{4,5}, Malachi Griffith, PhD^{1,2,3,6} and Obi L. Griffith, PhD^{1,2,3,6}

Purpose: Following automated variant calling, manual review of aligned read sequences is required to identify a high-quality list of somatic variants. Despite widespread use in analyzing sequence data, methods to standardize manual review have not been described, resulting in high inter- and intralab variability.

Methods: This manual review standard operating procedure (SOP) consists of methods to annotate variants with four different calls and 19 tags. The calls indicate a reviewer's confidence in each variant and the tags indicate commonly observed sequencing patterns and artifacts that inform the manual review call. Four individuals were asked to classify variants prior to, and after, reading the SOP and accuracy was assessed by comparing reviewer calls with orthogonal validation sequencing.

Results: After reading the SOP, average accuracy in somatic variant identification increased by 16.7% (p value = 0.0298) and average interreviewer agreement increased by 12.7% (p value < 0.001). Manual review conducted after reading the SOP did not significantly increase reviewer time.

Conclusion: This SOP supports and enhances manual somatic variant detection by improving reviewer accuracy while reducing the interreviewer variability for variant calling and annotation.

Genetics in Medicine (2018) <https://doi.org/10.1038/s41436-018-0278-z>

Keywords: somatic variant refinement; manual review

We are on a Coffee Break & Networking Session

Workshop Sponsors:

compute | calcul
canada | canada



Canadian Centre
for
Computational
Genomics

MicCM McGill initiative in
Computational Medicine